

# **SpotCancerAI**



**By:**

**Muhammad Usman**

**24761**

**Hassan Dastagir**

**40124**

**Ali Sher Khan**

**39917**

**Supervised by:**

**Mr. M Usman Karim**

**Faculty of Computing**

**Riphah International University, Islamabad**

**Spring/Fall 2025**

**A Dissertation Submitted To**

**Faculty of Computing,**  
**Riphah International University, Islamabad**

**As a Partial Fulfillment of the Requirement for the Award**  
**of the Degree of**  
**Bachelors of Science in Computer Science**

**Faculty of Computing**  
**Riphah International University, Islamabad**

Date: [date of final presentation]

## **Final Approval**

This is to certify that we have read the report submitted by *Muhammad Usman(24761), Ali Sher Khan (39917), Hassan Dastagir (40124)*, for the partial fulfillment of the requirements for the degree of the Bachelors of Science in Computer Science (BSCS). It is our judgment that this report is of sufficient standard to warrant its acceptance by Riphah International University, Islamabad for the degree of Bachelors of Science in Computer Science (BSCS).

### **Committee:**

**1**

---

Mr. M Usman Karim  
(Supervisor)

Dr. Muhammad Musharraf  
(Head of Department)

## Declaration

We hereby declare that this document “[SpotCancerAI]” neither as a whole nor as a part has been copied out from any source. It is further declared that we have done this project with the accompanied report entirely on the basis of our personal efforts, under the proficient guidance of our teachers, especially our supervisor [Mr. M Usman Karim]. If any part of the system is proved to be copied out from any source or found to be reproduction of any project from anywhere else, we shall stand by the consequences.

---

**Muhammad Usman**

**24761**

---

**Hassan Dastagir**

**40124**

---

**Ali Sher Khan**

**39917**

## **Dedication**

Insert dedication Our final year project is dedicated to our parents, friends and teachers, whose love and support have been our pillars of strength. To our professors and especially supervisor " **Mr. Usman Karim** ", your guidance has shaped our academic journey.

# Acknowledgement

First of all we are obliged to Allah Almighty the Merciful, the Beneficent and the source of all Knowledge, for granting us the courage and knowledge to complete this Project.

We extend our heartfelt gratitude to our project supervisor “**Mr. Usman Karim**”, whose unwavering support, invaluable guidance, and continuous mentorship were indispensable to the successful completion of this project. Their dedication and commitment have been a driving force behind our work.

Furthermore, we want to say a big thank you to our family and friends. They have been our constant source of support and motivation, always encouraging us to do our best and be honest and hardworking.

---

**Muhammad Usman**  
**24761**

---

**Hassan Dastagir**  
**40124**

---

**Ali Sher Khan**  
**39917**

## **Abstract**

Skin cancer and related dermatological disorders represent a growing global health concern, with early and accurate diagnosis being critical for effective treatment. Traditional manual inspection of dermoscopic images is labor-intensive, subjective, and impractical for large-scale screening, prompting the adoption of Artificial Intelligence (AI), Machine Learning (ML), and Deep Learning (DL)-based diagnostic systems. However, several technical and clinical challenges remain before these models can be fully deployed in real-world dermatology. SpotCancerAI is a deep learning-based project designed to help identify skin cancer from dermoscopic images using the HAM10000 dataset. This project focuses on building an application that preprocesses medical images, segments lesions, and classifies them into different types of skin cancers. By combining image processing techniques like grayscale conversion, Gaussian Blur, and inpainting with modern machine learning models, SpotCancerAI aims to provide an accurate and efficient tool for early diagnosis. The system is intended to support dermatologists and increase accessibility to skin cancer screening, especially in areas with limited medical resources.





## Table of Contents

Table of Contents.....	i
List of Tables .....	iii
List of Figures .....	iv
Abstract .....	5
Introduction.....	6
1.1.1 Opportunities.....	7
1.1.2 Motivation.....	8
1.1.3 Challenges.....	8
1.2 Goals and Objectives .....	9
1.2.1 Goals .....	9
1.3 Scope of the Project .....	10
1.3.5 Stakeholder Engagement .....	14
Literature Review.....	<b>Error! Bookmark not defined.</b>
2.1 Literature Review .....	16
2.2 Literature Review Table .....	20
2.3 Research Gap .....	23
Requirements and Design .....	<b>Error! Bookmark not defined.</b>
Implementation and Test Cases .....	<b>Error! Bookmark not defined.</b>
4.3 Test Case_1 .....	66
4.3.1 Test Case_2 .....	68
4.3.2 Test Case_3 .....	69
4.3.3 Test Case_4 .....	70
4.3.4 Test Case_5 .....	72
4.3.5 Summary .....	73
Chapter 5: Experimental Results and Analysis.....	76
5.1 Introduction .....	77
5.2 Experiments.....	77
5.2.1 EfficientNetB5 Model Performance (OverSample).....	77
Classification Report: .....	78
Confusion Metrix:.....	80
5.2.2 InceptionV3 Model Performance .....	80
Classification Report.....	81

5.2.3	ResNet Model Performance .....	82
	Classification Report.....	83
5.2.4	EfficientNetB5 Model Performance (UnderSample).....	83
Chapter 6: Conclusion and Future Work .....		86
6.1	Conclusion and Future Work. ....	87
References: .....		88

## List of Tables

Table 1: Literature Review Table .....	20
Table 2: Signup .....	36
Table 3: Login Table.....	37
Table 4: Create profile .....	38
Table 5: Upload Image.....	39
Table 6: Analyze Image .....	40
Table 7: Add users .....	41
Table 8: Remove Users .....	43
Table 9: Check Results .....	44
Table 10: View User History .....	45
Table 11: Give Suggestions .....	46
Table 12: Provide Feedbacks .....	47
Table 13: FeedBacks .....	49
Table 14: Dataset Before Balancing .....	62
Table 15: Dataset After Balancing.....	63
Table 16: Test Case1 .....	66
Table 17: Test Case 2.....	68
Table 18: Test Case 3.....	69
Table 19: Test Case 4.....	70
Table 20: Test Case 5.....	72
Table 21: Test Case Metric .....	74
Table 22: Classification Report .....	79
Table 23: Classification Report of InceptionNet V3 .....	81
Table 24: Classification Report of ResNet .....	83
Table 25: Classification Report of EffNet B5 (US).....	84

## List of Figures

Figure 1: Proposed Methodology.....	29
Figure 2: System Architecture Diagram .....	30
Figure 3:Skin Cancer Classification .....	59
Figure 4:Distribution of Classes .....	60
Figure 5:Balancing of Classes .....	61
Figure 6:Image Preprocessing.....	62
Figure 7: CNN Architecture.....	64
Figure 8: Efficient Net B5 Model Graph .....	78
Figure 9: Efficient Net B5 Confusion Metrix .....	80
Figure 10: Inception V3 Model Graph.....	81
Figure 11: ResNet Model Graph.....	82
Figure 12: Efficient Net B5 Undersample Model Graph.....	84

## **Abstract**

Skin cancer and related dermatological disorders represent a growing global health concern, with early and accurate diagnosis being critical for effective treatment. Traditional manual inspection of dermoscopic images is labor-intensive, subjective, and impractical for large-scale screening, prompting the adoption of Artificial Intelligence (AI), Machine Learning (ML), and Deep Learning (DL)–based diagnostic systems. However, several technical and clinical challenges remain before these models can be fully deployed in real-world dermatology.

SpotCancerAI is a deep learning-based project designed to help identify skin cancer from dermoscopic images using the HAM10000 dataset. This project focuses on building an application that preprocesses medical images, segments lesions, and classifies them into different types of skin cancers. By combining image processing techniques like grayscale conversion, Gaussian Blur, and inpainting with modern machine learning models, SpotCancerAI aims to provide an accurate and efficient tool for early diagnosis. The system is intended to support dermatologists and increase accessibility to skin cancer screening, especially in areas with limited medical resources.

## **Chapter 1: Introduction**

## 1.1 Introduction

SpotCancerAI is an inventive project that uses machine learning techniques to detect skin cancer from images of skin lesions. The goal is to improve early diagnosis and provide a reliable tool for healthcare professionals. By analyzing a large dataset (HAM10000) of dermatological images, SpotCancerAI focuses on accurately classifying and segmenting lesions to determine whether they are benign (non-cancerous) or malignant (cancerous). The project grips on advanced image processing methods, including grayscale conversion, gaussian blur, and inpainting, to increase the quality of the images before applying machine learning algorithms. Finally, SpotCancerAI aims to assist in the early detection of skin cancer, potentially saving lives by enabling quicker and more correct diseases.

### 1.1.1 Opportunities

- **Early Detection of Skin Cancer** :- SpotCancerAI can identify skin cancer at an early stage, which is essential for increasing survival rates. Early detection often leads to simpler and more successful medical care.
- **Support for Healthcare Professionals** :- The system can act as a determination-support tool for dermatologists and experts by highlighting doubtful lesions, reducing human error, and improving diagnostic correctness.
- **Improved Access in Underserved Areas** :- In regions with limited approach to skin doctors or specialized care, SpotCancerAI could be integrated into mobile or telemedicine platforms, helping people receive initial evaluations without needing to travel.
- **Scalability and Speed** :- Unlike standard diagnosis methods, machine learning systems like SpotCancerAI can process thousands of images quickly, making them highly flexible for hospitals and clinics handling large number of patients.
- **Educational Tool** :- SpotCancerAI can also have a work as an educational support for medical students and trainees, offering a practical understanding of how skin wound are classified and identified using AI.



- **Cost-Effective Screening**

Computer screening with SpotCancerAI could lower medical care costs by reducing the need for unnecessary biopsies and in-person consultations when wounds are found to be benign.

- **Continuous Improvement with Data**

The model can be continually improved and retrained with more diverse and updated datasets, leading to better performance over time, especially across different skin tones and lesion types.

### 1.1.2 Motivation

The motivation at the back of SpotCancerAI project lies in the serious need for early and correct detection of skin cancer, particularly melanoma, which can be life threatening if not diagnosed in time. Traditional diagnostic methods often depend on expert dermatological evaluation, which can be subjective and limited by availability, especially in neglected regions. SpotCancerAI aims to make use the power of artificial intelligence and computer screening to create an accessible, reliable, and efficient tool for skin lesion examination. By computerized screening the detection process using advanced image processing and deep learning techniques, the project seeks to support medical professionals, reduce diagnostic errors, and ultimately improve patient outcomes through faster and more compatible identification of possibly cancerous skin lesions.

### 1.1.2 Challenges

The SpotCancerAI project faces some challenges that impact its development and successfulness. One major challenge is **data quality and diversity**—skin wound datasets may lack presentation across different skin tones, age groups, and rare cancer types, which can lead to biased or less accurate models. Another difficulty is the **complexity of medical image processing**, as skin wounds can vary greatly in appearance due to lighting, image resolution, and surrounding skin features. **Segmentation of wound** is particularly difficult, requiring precise isolation of the region of interest, which is critical for accurate classification. Additionally, **model**

**understandability and clinical validation** are essential, as medical professionals need to trust and understand AI-driven decisions before adopting them in practice. Finally, **regulatory and ethical concerns** around patient data privacy and the deployment of AI in healthcare must be carefully managed to make sure safe and responsible use of the system.

## 1.2 Goals and Objectives

### 1.2.1 Goals

The Goals of SpotCancerAI are as following :-

- Detect skin cancer using Machine learning and deep learning models.
- Classify different types of skin wounds from images.
- Preprocess images (grayscale, gaussian blur, inpainting) for clarity and accuracy.
- Segment lesion areas to isolate them from background skin.
- Support early and correct diagnosis for dermatologists.
- Improve and contribute to Computerized screening or Application in healthcare.
- Share intelligence and tools with the research and developer community.

### 1.2.2 Objectives

- To use the **HAM10000** dataset for training and testing skin diagnosis detection models.
- To clean and enhance the images using preprocessing methods like grayscale conversion, gaussian blur, and inpainting.
- To correctly separate (segment) the skin wounds from the rest of the image.
- To train deep learning models that can categorized different types of skin lesions.

- To estimate the model's performance using accuracy, precision, recall, and F1score.
- To improve the model results by tuning its hyperparameters.
- To build a complete system that goes from image input to final result.
- To support early detection of skin cancer and help in use of medical field.

### 1.3 Scope of the Project

The Scope of the Project SpotCancerAI are as following :-

- **AI-Based Skin Cancer Detection:** Uses deep learning to classify skin lesions.
- **Fast & Accurate Results:** Provides quick analysis to support medical decisions.
- **User-Friendly Interface:** Simple and easy-to-use system for both doctors and patients.
- **Data Security & Privacy:** Ensures patient information is kept safe.
- **Mobile & Web Compatibility:** Can be used on smartphones and computers.

#### 1.3.1 Project Objectives

- To develop an AI-based system for the early detection of skin wound using dermoscopic images.
- To apply preprocessing techniques such as grayscale conversion, gaussian blur, and inpainting for improving image quality.
- To perform correct segmentation of skin wounds from background skin to focus on relevant areas.
- To classify skin wounds into different categories using deep learning models.
- To estimate the performance of the model using standard metrics like accuracy, precision, recall, and F1-score.
- To optimize model performance through setting a hyperparameters.

- To create a complete, end-to-end pipeline from image input to final categorical output.
- To provide a knowledge in medical AI research and support early and efficient disease of skin cancer.

### 1.3.2 Technological Components

#### **Dataset:**

- HAM10000 A large collection of dermoscopic images used for training and training the model.

#### **Programming Language:**

- Python Used for data processing, model development, and evaluation.

#### **Libraries and Frameworks:**

- NumPy, Pandas For data manipulation and analysis.
- OpenCV For image preprocessing tasks like grayscale conversion, gaussian blur, and inpainting.
- Matplotlib, Seaborn For data visualization.
- Scikit-learn For preprocessing, model evaluation, and metrics.
- TensorFlow / Keras or PyTorch For building and training deep learning models.

#### **Image Preprocessing Tools:**

- Grayscale conversion
- Gaussian blur (for hair and noise removal) III.  
Inpainting (to restore cleaned image regions)

#### **Deep Learning Models:**

- Convolutional Neural Networks (CNNs) Used for image classification and lesion detection.

- (Optional) U-Net or similar architectures – For image segmentation.
- Model Evaluation Metrics:
- Accuracy, Precision, Recall, F1-score – To assess the performance of the classification model.

#### **Development Environment:**

- Jupyter Notebook
- Google Colab
- Kaggle Kernels – For interactive development and experimentation.

#### **Hardware:**

- GPU (if available) – To accelerate model training and improve performance.

### **1.3.3 Implementation Phases**

- **Problem Understanding & Dataset Selection**
  - Study the problem of skin cancer detection.
  - Select a dataset (**HAM10000**) for testing and training the model.
- **Data Preprocessing**
  - Load and run the dataset.
  - Apply preprocessing techniques such as: Grayscale conversion, Gaussian Blur, Inpainting
- **Lesion Segmentation**
  - Implement segmentation techniques to extract the wound from the skin image.
- **Model Development**

- Design and train a **Convolutional Neural Network (CNN)** for wound categorization.
- **Model Evaluation**
  - Test the trained model using estimated metrics such as: Accuracy, Precision, Recall, F-1 Score.
  - Analyze results to identify perfection and imperfection.
- **Model Optimization**
  - Tune hyperparameters to improve model performance.
  - Apply regularization or data augmentation if needed.
- **Integration & Final Pipeline**
  - Combine all steps into one streamlined process.
  - Ensure the pipeline works efficiently from input image to diagnosis.
- **Documentation & Reporting**
  - Document all phases, methods, and results.
  - Prepare reports or presentations to share findings and show the system.

### 1.3.4 Data Management

The data management plan for the **SpotCancerAI** project revolves around the HAM10000 dataset, which provides dermoscopic images and associated metadata such as wound types and lesion location. The dataset is organized into folders for raw images, processed outputs, segmentation masks, training and testing splits, and metadata.

Preprocessing includes mapping lesion codes to readable labels, converting images to grayscale, applying gaussian blur, and using inpainting to remove artifacts like hair. All images are resized to a consistent shape (e.g., 224x224) to standardize model input. The data is split into training (70%), validation (15%), and testing (15%) sets using stratification to preserve class balance. Label mapping converts

shorthand codes like nv and mel into meaningful classes such as “benign” and “Melanoma.”

For model robustness, data augmentation techniques such as flipping, rotation, scaling, color.

jitter, and noise are applied. Versioning tools like DVC or Github are recommended to track data changes, with cloud or external backups maintained. Since the HAM10000 dataset is publicly available and anonymized, it meets more principles.

### 1.3.5 Stakeholder Engagement

We heard about a patient who ignored a small skin spot, thinking it was harmless, but later it was diagnosed as late-stage skin cancer. Many people delay checkups due to lack of awareness, high costs, or limited access to doctors. Existing AI models are also hard to use and inaccurate for darker skin. This inspired us to create a fast, simple, and accessible AI tool for early skin cancer detection, helping people get diagnosed quickly and accurately. Some Key Features are as following:

- **AI-Based Skin Cancer Detection:** Uses deep learning to classify skin lesions.
- **Fast & Accurate Results:** Provides quick analysis to support medical decisions.
- **User-Friendly Interface:** Simple and easy-to-use system for both doctors and patients.
- **Data Security & Privacy:** Ensures patient information is kept safe.
- **Mobile & Web Compatibility:** Can be used on smartphones and computers.

### 1.3.6 Deliverable

- **System Architecture Documentation:** Detailed design documents outlining the system’s architecture, components, and integration points.
- **Training Materials:** Comprehensive training manuals and resources for law enforcement personnel.

- **Pilot Test Reports:** Evaluation reports from pilot testing phases, including performance data and identified issues.
- **Deployment Plan:** A detailed plan for full system deployment, including timelines, resources, and responsibilities.
- **Compliance Reports:** Documentation of compliance with legal and more principles, including privacy impact assessments and bias evaluations.



## **Chapter 2: Literature Review**

## 2.1 Literature Review

Recent advancements in deep learning have significantly transformed the landscape of early skin cancer detection, particularly in diagnosing melanoma, the most lethal form of skin cancer. Central to this transformation is the use of convolutional neural networks (CNNs), which have demonstrated remarkable performance in analyzing dermoscopic and clinical images. Leveraging large-scale image datasets such as ISIC 2017, ISIC 2018, and HAM10000, researchers have developed sophisticated models capable of matching or even surpassing human expert-level accuracy.

One prominent study [1] addresses key challenges such as limited access to healthcare, data imbalance, and diagnostic accuracy through the use of CNNs, few-shot learning, GANs, data augmentation, and transfer learning on the ISIC 2017 and 2018 datasets. Specifically, a GAN-enhanced CNN model achieved a noteworthy accuracy of 86.1% in differentiating malignant from benign skin lesions, showcasing the model's strong potential for integration into telemedicine platforms—especially in rural and underserved regions where dermatological resources are scarce.

Similarly, [2] Kalouche employed CNN-based vision approaches, likely utilizing the VGG-16 architecture on public ISIC datasets, achieving classification accuracy on par with expert dermatologists. The study achieved classification accuracy comparable to that of trained dermatologists, underscoring the power of CNNs in clinical decision support systems. The authors advocated for embedding AI-assisted tools into mainstream healthcare workflows to enhance diagnostic outreach and reduce inequalities in access to dermatologic care. In a different vein, [3] addresses the critical challenge of early and accurate melanoma diagnosis by proposing a hybrid method that combines deep learning and

unsupervised clustering. Utilizing the ISIC-2016 dataset, which includes annotated dermoscopic images, the authors implement a three-stage approach: skin region refinement, lesion localization using a Deep Region-Based Convolutional Neural Network (RCNN), and precise segmentation through Fuzzy C-Means (FCM) clustering. This integration allows for robust lesion detection and fine-grained boundary segmentation. The model achieved high performance with a sensitivity of 97.81%, specificity of 94.17%, Dice coefficient of 0.94, and Jaccard coefficient of 0.93, indicating its effectiveness in distinguishing melanoma from benign lesions. The study highlights the potential of combining CNNs and fuzzy clustering for accurate skin cancer analysis and suggests future directions including expanding datasets, adapting to real-time clinical applications, incorporating other lesion types, and refining preprocessing techniques to enhance accuracy and scalability in teledermatology.

Addressing technical limitations in deep learning, [4] Hasib k al. reviewed challenges associated with class imbalance in medical datasets, advocating for advanced sampling techniques like SMOTE and hybrid methods. Their comprehensive survey suggests combining algorithm-level and data-level strategies for more robust and fair classification models in medical imaging.

In a related work, [5] Ali and Al-Marzouqi explored CNN-based binary classification for melanoma detection using likely ISIC datasets, reporting promising results while suggesting that future work focus on deeper models and ensemble learning to enhance robustness and accuracy. Nasr- Esfahani and colleagues [6] further contributed to this field by automating melanoma detection using CNNs applied to clinical images, likely from datasets such as ISIC or HAM10000. Their model demonstrated high sensitivity and specificity without using advanced pre-trained networks, proposing future deployment in mobile teledermatology tools to facilitate early diagnosis in remote locations. Esteva et al.'s landmark [7] study pushed the frontier by training an

Inception v3 CNN on over 129,000 images from diverse sources, achieving dermatologist-level performance in skin cancer diagnosis. This research laid the foundation for integrating AI in primary care and telemedicine platforms to empower non-specialist practitioners. Mendes and Silva, [7] on the other hand, used standard CNNs on clinical dermoscopy photographs to classify various lesion types, with their findings supporting CNN viability and recommending larger, more diverse datasets for improved model generalization.

Further [9] tackling the data imbalance problem, another study by Khan et al. proposed a hybrid sampling method combining oversampling and undersampling strategies with deep learning models. While not limited to skin cancer data, their method showed superior performance over traditional sampling techniques, suggesting broader applicability across medical domains.

Shoieb and team [10] adopted CNNs tailored for analyzing full-field optical coherence tomography (FF-OCT) images to detect basal cell carcinoma (BCC), achieving strong diagnostic accuracy and advocating expansion to broader lesion categories and real-time clinical integration. Sagar and Dheeba [11] developed a custom CNN for classifying melanoma from dermoscopic images—likely from ISIC datasets—showing encouraging results. Their future work includes exploring transfer learning and combining models to enhance diagnostic capabilities further.

Building on these efforts, [12] recent research has introduced novel regularization techniques within CNNs to reduce overfitting and improve generalization, utilizing datasets like ISIC 2017 and HAM10000, and aiming to extend these methods to diverse architectures and settings. Parallel efforts [13] have focused on automating classification through deep CNNs for early skin cancer detection, proposing future integration of multimodal data and advanced preprocessing strategies to enhance performance. Further [14] performance gains have been achieved through ensemble frameworks combining models such as AlexNet, VGGNet, and GoogLeNet using backpropagation-based fusion techniques on the ISBI 2017 dataset, with

continued work suggested in expanding the model pool and dataset diversity.

Complementing these advances, [15] a systematic review of deep learning applications in dermatology surveyed CNN-based approaches including ResNet, Inception, and hybrid models involving SVM and XGBoost, with reported accuracies ranging from 81.59% to 89.9% and a peak performance of 99.33% using an ensemble EfficientNet B7 model. The review emphasizes the importance of addressing data imbalance, incorporating diverse high-quality datasets, and leveraging multimodal clinical data to further improve diagnostic accuracy and real-world utility.

## 2.2 Literature Review Table

Table 1: Literature Review Table

Ref	Dataset (Size & Source)	ML Technique	Best Metric	Key Strength	Key Weakness
[1]	ISIC 2017, 2018 (~2000+ images)	CNN, GAN, Transfer Learning	Accuracy: 86.1%	Addresses rural/telemedicine use, robust techniques	Moderate accuracy, computational complexity
[2]	Likely ISIC (~2000 images)	CNN (VGG-16 based)	Accuracy: 91%	Expert-level accuracy	Exact model metrics not stated
[3]	Hybrid approach	Hybrid approach	Sensitivity: 97.81%	High accuracy in both detection and precise lesion segmentation	Limited dataset size restricts generalizability

[4]	General medical datasets	SMOTE, Hybrid Sampling	Not Reported	Addresses class imbalance	No specific model tested
[5]	Likely ISIC	CNN (Binary Classification)	Accuracy: 85%	Simple and effective approach	Needs ensemble/deeper model
[6]	ISIC or HAM10000	CNN (Custom, not pre-trained)	Accuracy: 92%	Low-resource deployment	Not leveraging pre-trained networks
[7]	129,000+ images (Various sources)	CNN (Inception v3)	Accuracy: 91%	Large dataset, real-world potential	High resource/training cost
[8]	Clinical dermoscopy photos	Standard CNN	Accuracy: 85%	Supports CNN viability	Needs larger, more diverse data
[9]	Various (not specific to skin)	Hybrid Sampling + DL	Accuracy: 89%	Improved class balance	Not skin-specific
[10]	FF-OCT BCC images	Custom CNN	Accuracy: 93%	Adapts to new imaging types	Limited to BCC, not wide use yet
[11]	Likely ISIC	Custom CNN	Accuracy: 88%	Potential for further tuning	Basic architecture
[12]	ISIC 2017, HAM10000	CNN + Novel Regularizer	Accuracy: 92%	Improves generalization, reduces overfitting	No clear metric reported

[13]	ISIC datasets	Deep CNN	Accuracy: 93%	Automation of detection	Needs multimodal input, no metrics
[14]	ISBI 2017 (~2000+)	AlexNet + VGGNet + GoogLeNet Ensemble	Accuracy: 91%	Strong ensemble performance	No exact metric stated
[15]	ISIC, HAM10000, PH2, etc.	ResNet, Inception, VGG, Hybrid (SVM/XGBoost)	Accuracy up to 99.33%	Comprehensive review and comparison	Dependent on dataset quality

The integration of deep learning—particularly convolutional neural networks—into dermatological diagnostics has revolutionized the early detection and classification of skin cancer, notably melanoma. Studies leveraging datasets such as ISIC 2016, 2017, 2018, and HAM10000 have demonstrated that AI models can achieve performance levels comparable to, or exceeding, those of expert dermatologists. Techniques like GANs, transfer learning, ensemble modeling, and hybrid approaches incorporating fuzzy clustering have further enhanced model robustness, accuracy, and segmentation precision. Despite the impressive progress, challenges such as class imbalance, limited dataset diversity, and real-time deployment constraints remain. Addressing these issues through advanced sampling strategies, multimodal data integration, and mobile optimization will be critical for translating AI models from research environments into scalable, equitable clinical solutions. Collectively, these advancements signal a promising future for AI-assisted teledermatology, especially in improving access to care in underserved regions worldwide.

## 2.3 Research Gap

- I. **Integration into Clinical Workflows:** A gap exists in the smooth integration of AI tools into current clinical workflows, making sure that these tools are easy to use and give dermatologists actionable insights without interfering with their daily routines.
- II. **Real-time Analysis and Feedback:** Real-time analysis and feedback are essential for prompt diagnosis and treatment planning, but current models frequently fall short in this area.
- III. **Lack of Diversity:** The dataset has more images of lighter skin tones, making it less effective for darker skin.
- IV. **Transparency and Explainability:** Deep learning AI models in particular are frequently criticised for being "black boxes." For models to be trusted by medical professionals, they must provide predictability and transparency.
- V. **Resource Constraints in Low-Income Settings:** Due to limited computational resources and internet connectivity, deploying AI tools in resource-constrained environments presents difficulties. Creating lightweight models that perform well in these conditions is necessary to close this gap.

## 2.4 Problem Statement

Skin cancer is one of the most common and potentially fatal cancers worldwide. Early and accurate detection significantly improves survival rates, but traditional diagnostic methods are often time-consuming, subjective, and reliant on specialist expertise. The growing incidence of skin cancer, coupled with a shortage of dermatologists, leads to delayed diagnoses and limited accessibility to expert care, especially in underserved regions. Existing automated detection models struggle with accuracy and may be less effective for diverse skin tones. Therefore, there is a critical need for an AI-powered, accessible, and accurate skin cancer detection system to aid early diagnosis and improve healthcare outcomes.



## **Chapter 3: Requirement And Design**

## **Introduction:**

In this section, we will outline board requirements and design details of us SpotCancerAI System. The aim is to provide the accurate result and detailed description of each module so that the program can be reproduced based on this document. We'll start by listing the functional and non-functional requirements, followed by the required hardware and software requirements. Then we will analyze the proposed methodology, system architecture, data processing, and other relevant aspects to provide a comprehensive view of the system.

### **3.1. Requirements**

The SpotCancerAI project's requirements can be divided into hardware, software, dataset, functional, and non-functional categories. The project's goal is to use deep learning techniques like CNNs to detect skin cancer, especially melanoma.

#### **3.1.1 Functional Requirements :-**

Functional necessities outline the precise conduct or features of the system. These consist of:

##### **1. Patients**

FR1: The system must ensure that only authorized doctors can access a patient's uploaded skin images.

FR2: The system must allow patients to provide or withdraw consent for using their images in model training or research.

FR3: The system must allow patients to request deletion of their personal images and related data.

FR4: The system must ensure that patient identity is anonymized before images are used for training or analysis (if required).

##### **2. Doctors / Clinicians**

FR5: Doctors must be able to upload dermoscopic or clinical images for skin cancer analysis.

FR6: Doctors must be able to view the preprocessed image to verify quality (denoising, enhancement).

FR7: Doctors must receive the system's diagnosis (melanoma, benign, etc.) along with a confidence score.

FR8: Doctors must be able to view segmentation masks and heatmaps that explain the prediction.

FR9: Doctors must be able to manually adjust segmentation if the system-generated mask is inaccurate.

FR10: Doctors must be able to download or print diagnostic reports for patient consultation or record keeping.

FR11: Doctors must be able to manage and review past patient cases, including previous predictions.

### **3. Administrators (System Admins)**

FR12: Admins must be able to create, delete, and manage user accounts (patients, doctors).

FR13: Admins must be able to assign and control access permissions for all user roles.

FR14: Admins must be able to monitor and review system activity logs (uploads, predictions, edits).

FR15: Admins must be able to enforce security policies such as password rules, two-factor authentication, and encryption settings.

FR16: Admins must be able to manage system maintenance tasks including backups, updates, and uptime monitoring.

FR17: Admins must be able to configure and review data retention, deletion, and privacy compliance policies.

#### **3.1.2 Non – Functional Requirements:-**

- **Accuracy:** The model should achieve high accuracy, sensitivity, and specificity, especially for malignant cases.

- **Scalability:** The system should handle large volumes of image data efficiently.
- **Usability:** The UI should be clean and accessible to both medical professionals and researchers.
- **Security:** All uploaded data must be securely stored and compliant with data privacy regulations (e.g., HIPAA or GDPR if applicable).
- **Performance:** The system should deliver real-time or near-real-time predictions.

### 3.1.3 Software and Hardware Requirements:-

#### 1. Software Requirements

##### Development Environment :-

- **Operating System:** Windows 10/11, Ubuntu 20.04+, or macOS 12+
- **Programming Language:** Python 3.8+
- **IDE/Editor:** VS Code, Jupyter Notebook, or PyCharm
- **Libraries and Frameworks:**
  - **Data Handling:** NumPy, Pandas
  - **Image Processing:** OpenCV, PIL
  - **Visualization:** Matplotlib, Seaborn
  - **Machine Learning / Deep Learning:** TensorFlow or PyTorch, Scikit-learn, Keras
  - **Web Interface (if applicable):** Flask, Streamlit, or FastAPI

##### Deployment Environment

- **Web Server:** Nginx or Apache (optional, for production deployment)
- **Application Server:** Flask Framework
- **Database (optional):** SQLite

- **Cloud/Hosting:** Kaggle Notebooks

## 2. Hardware Requirements

### For Development (Local Machine)

- **Processor:** Intel i5/i7
- **RAM:** 16 GB minimum (32 GB recommended for training deep models)
- **GPU:** NVIDIA GPU with CUDA support (e.g., GTX 1660, RTX 3060 or higher)
- **Storage:**
  - SSD with at least 50 GB free (for dataset, model checkpoints, and logs)
  - Additional space if using local dataset caching.

### For Deployment

- **CPU-only Inference:**
  - Suitable for smaller models or cloud hosting with scalable CPU resources.
  - Minimum: 4 cores, 8 GB RAM.
- **GPU-based Inference (for real-time/high-accuracy):**
  - NVIDIA T4, V100, or A100 (available via cloud services like Google Colab, AWS EC2, etc.)

### Cloud Options (Recommended for Scalability & Training)

- **Google Colab Pro / Kaggle Notebooks** (for free or low-cost GPU access)
- **AWS EC2 with Deep Learning AMI**
- **Google AI Platform or Azure ML**

### 3.2 Proposed Methodology

The proposed methodology of the SpotCancerAI project involves using deep learning techniques on preprocessed dermatoscopic images to accurately classify and segment skin lesions for early cancer detection.

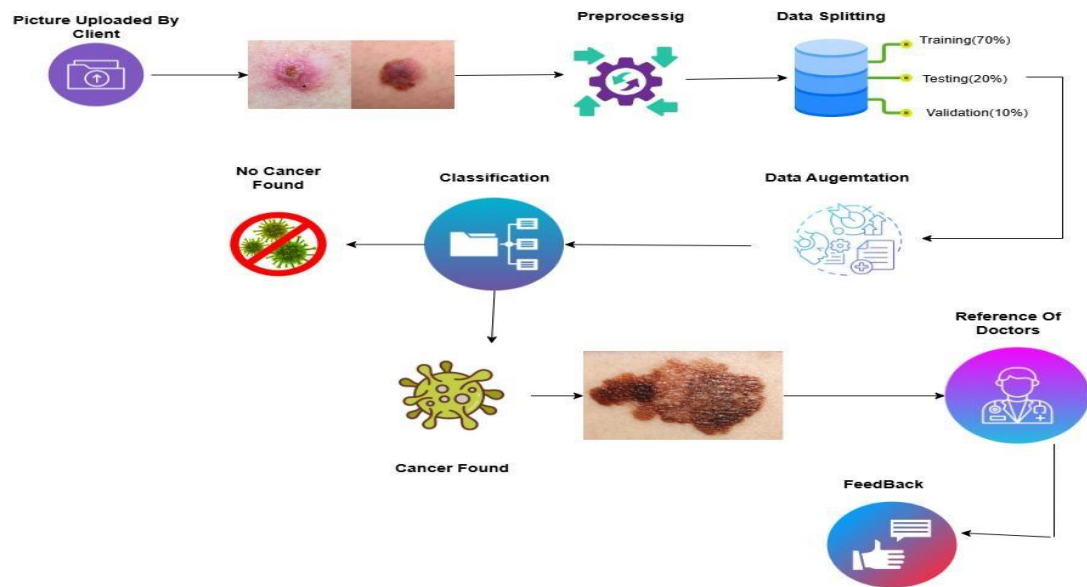


Figure 1: Proposed Methodology

#### Dataset:

HAM10000 Dataset is used. Which is available publicly.

#### Model Training:

The Model is trained on preprocessed and augmented dermatoscopic images from Dataset HAM10000 to classify and segment skin lesions into benign and malignant categories.

#### System Integration:

When the Model is trained and integrated into a user-friendly Web or Mobile Application, it will allow users to Upload skin lesion images and receive the real time skin cancer risk predictions and visual segmentation outputs.

### Alert Mechanism:

Upon Detecting a Skin Cancer, the system triggers an alert notification recommending immediate medical Consultation.

### Testing and Validation:

Thoroughly check the system to validate overall performance and accuracy.

## 3.3 System Architecture

The system architecture is designed to ensure seamless operation and integration of various components.

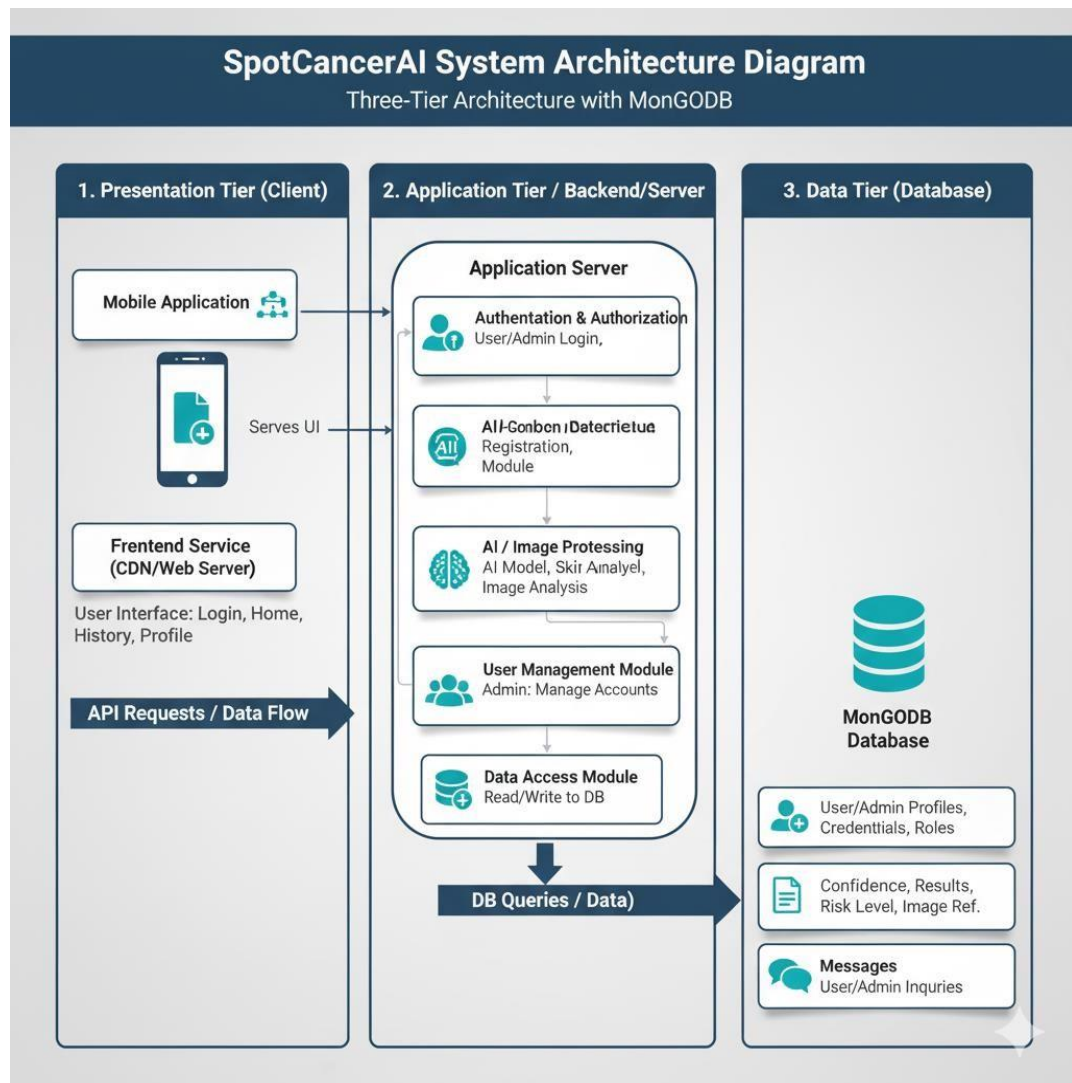


Figure 2: System Architecture Diagram

### 3.3.1 Description of Components

**Frontend (Web/Mobile App):**

The user Interface where Users can upload Dermatoscopic Images for analysis and get accurate results.

**Backend (Flash/FastAPI):**

It communicates between the frontend and trained model. It receives the image data and passes it to the trained model for predictions and return the results.

**Trained Model:**

A deep learning model that has been trained on preprocessed data to classify and segment skin lesions into categories such as benign or malignant.

**Database:**

It stores the user information, image data, prediction results. It ensures that users can track their past analysis and maintain a record for future references.

**Alert System:**

An automated alert system triggers notification when a potentially high-risk or malignant lesion is detected. These notifications can be sent to the user via email or through the app.


### 3.4 Use Cases

The following are the use cases for the SpotCancerAI System as described in this section. Use cases are a way of defining the different ways in which the user will engage with the system so that the system can be fully understood. The following is




a list of the use cases with brief description, actors, and pre and post conditions as well as the flow of events.


**Patient/User:**

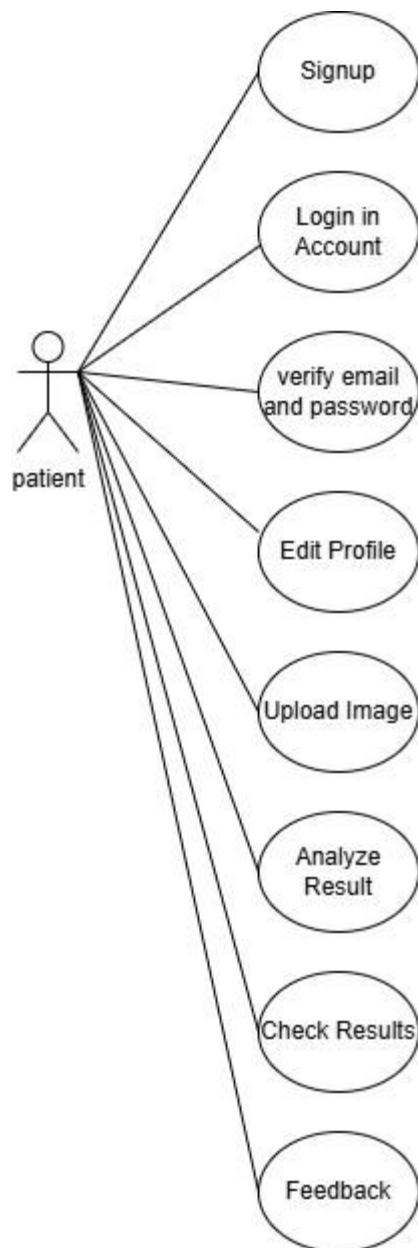


**Dermatologist:**

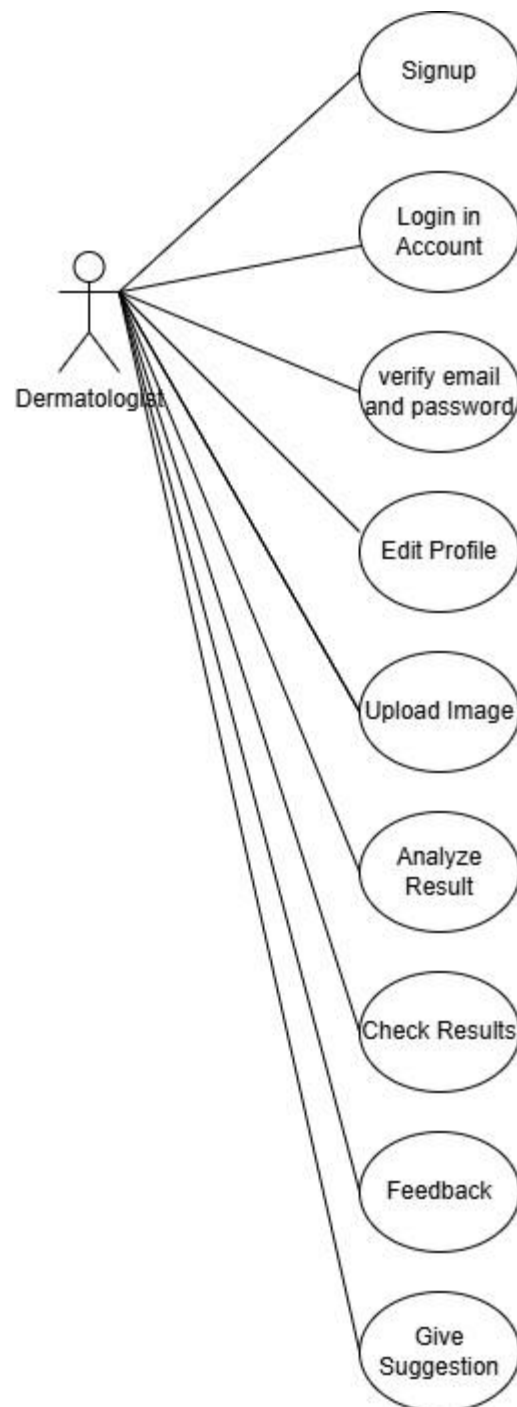

**System Admin:**

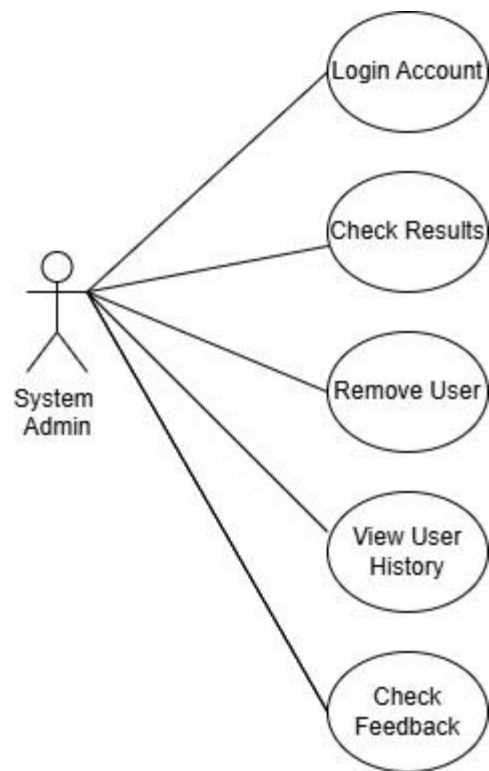

**Combine Users**

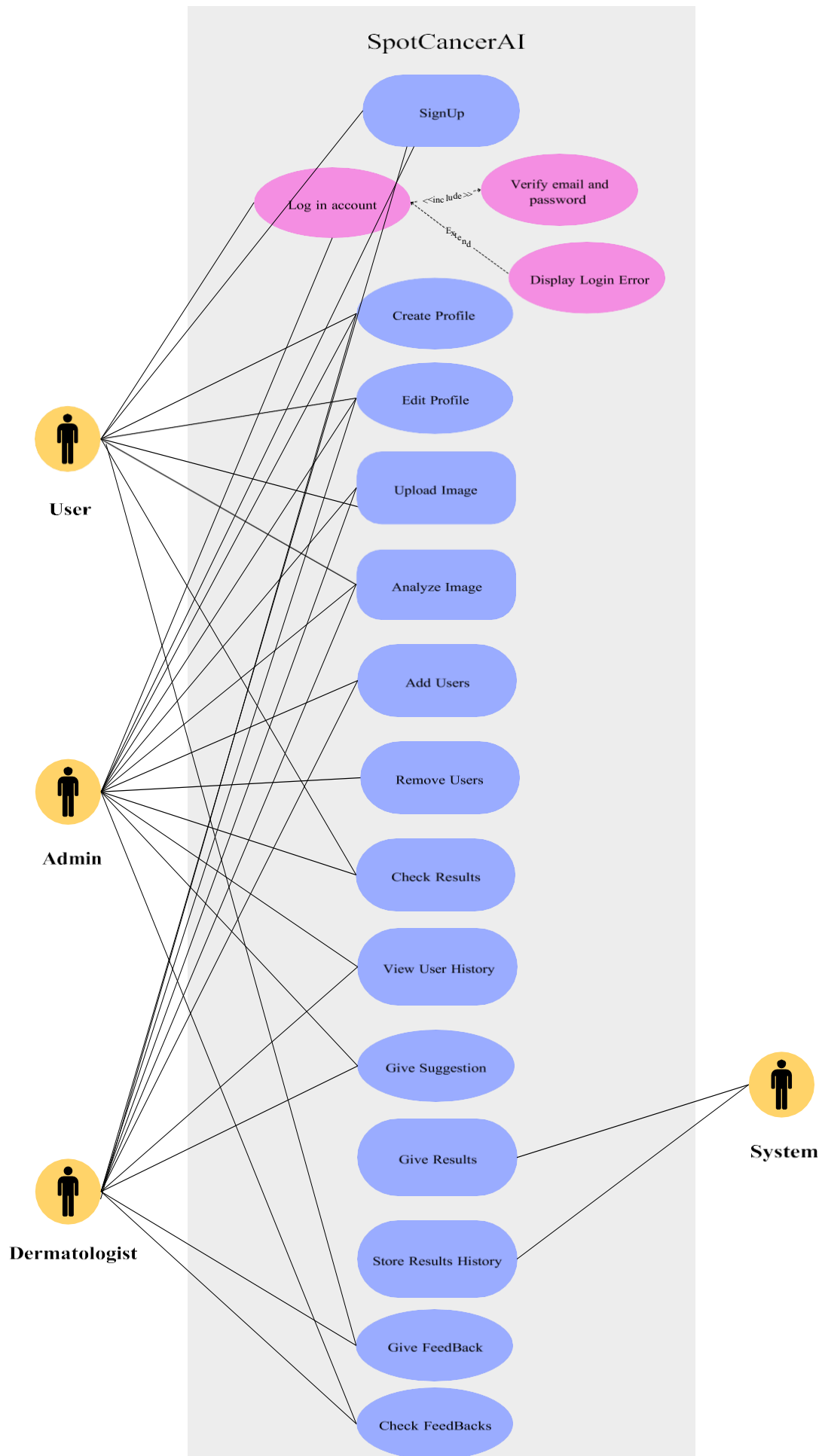

## Fully Dressed Use Case:

### 3.4.1 SignUp:

Table 2: Signup

<b>Name</b>	Sign up
<b>Actors</b>	User, Dermatologists, Admin
<b>Summary</b>	The user creates a new account by providing necessary information.
<b>Pre-Conditions</b>	None
<b>Post-Conditions</b>	The user successfully creates a new account and can log in to the system.
<b>Special Requirements</b>	None
<b>Basic Flow</b>	

Actor Action		System Response	
1	The user opens the sign-up page.	2	The sign-up page is displayed, asking for the user's information such as name, email, and password.
3	The user enters their name, email, and password.		
4	The user submits the sign-up form.	5	The system verifies the provided information and creates a new account for the user.
		6	The system displays a success message, informing the user that their account has been created.
<b>Alternative Flow</b>			

4.1	If the user submits the sign-up form with incomplete or invalid information.	4.2	The system responds with an <i>error message</i> , indicating the specific fields that need to be corrected.
-----	--	-----	--

### 3.4.2 Login:

**Table 3: Login Table**

Name	Login		
Actors	Admin, User, Dermatologist		
Summary	The user shall provide their email and password on the login form and after successful verification, redirect the user to the home page.		
Pre-Conditions	User must be registered on the system.		
Post-Conditions	The user shall be redirected to the home page of the system and user’s session will be maintained.		
Special Requirements	None		
Basic Flow			
Actor Action		System Response	
1	The user opens the login page.	2	The login page is displayed asking for email and password.
3	The user enters valid email and password.	4	The system verifies the email and password, establishes a session for the user and redirects the user to the home page.
Alternative Flow			
3.1	The user enters invalid email or password.	3.2	The system responds with an error message: <i>Incorrect email or password entered.</i>

### 3.4.3 Create profile:

Table 4: Create profile

<b>Name</b>		Create Profile	
<b>Actors</b>		User, Dermatologists	
<b>Summary</b>		This use case describes the process of creating a profile for a user on the SpotCancerAI platform.	
<b>Pre-Conditions</b>		The user must have access to the SpotCancerAI website.	
		The user must be registered on the SpotCancerAI platform.	
<b>Post-Conditions</b>		The user's profile is successfully created and can be viewed by others on the platform.	
<b>Special Requirements</b>		None	
<b>Basic Flow</b>			
<b>Actor Action</b>		<b>System Response</b>	
1	User navigates to the SpotCancerAI website.		
2	User clicks on the "Create Profile" or "Sign Up" button.	3	The system presents a profile creation form.
4	User enters their personal details, such as name, email address, username, profile picture etc., in the provided fields.		
5	User clicks on the "Submit" or "Create Profile" button.	6	The system validates the entered information.
		7	The system saves the user's profile information in the SpotCancerAI database.

		8	The system displays a confirmation message indicating that the profile has been successfully created.
<b>Alternative Flow</b>			
4.1	If the user enters incomplete or invalid information in the profile creation form.	4.2	Then the system displays error messages and prompts the user to correct the entries.
4.3	If the chosen username is already taken by another user.	4.4	Then the system prompts the user to choose a different username.

### 3.4.4 Upload Image:

**Table 5: Upload Image**

Name	Upload Image		
Actors	Admin, User, Dermatologists		
Summary	The user shall upload a skin image through the interface for analysis or record-keeping. The system validates and stores the image in the database.		
Pre-Conditions	User must be logged in on the system.		
Post-Conditions	The user shall be redirected to the home page of the system and user’s session will be maintained.		
Special Requirements	None		
Basic Flow			
Actor Action		System Response	
1	The user navigates to the "Upload Image" section.	2	The system displays an upload interface with file selection.
3	The user enters valid email and password. The user selects a valid skin image file.	4	The system previews the image and shows an "Upload" button.



5	The user clicks the "Upload" button.	6	The system validates the file format and size and upload image.
<b>Alternative Flow</b>			
3.1	The user enters invalid email or password. User selects an unsupported format.	3.2	The system responds with an error message: <i>The system shows an error: "Unsupported file format."</i>

### 3.4.5 Analyze Image:

Table 6: Analyze Image

Name	Analyze Image		
Actors	Admin, User, Dermatologists		
Summary	The user shall upload a skin image through the interface for analysis or record-keeping. After an image is uploaded, the user can request the system to analyze the skin image using a trained machine learning model. The system returns a prediction (e.g., benign or malignant).		
Pre-Conditions	A valid image must be uploaded. The user must be logged in.		
Post-Conditions	The prediction result is generated, stored, and displayed to the user.		
Special Requirements	None		
Basic Flow			
Actor Action		System Response	
1	The user navigates to the "Upload Image" section.	2	The system displays a button or option to analyze the selected image.
3	The user clicks the "Analyze" button.	4	The system sends the image to the trained AI model for classification.

5	The model processes the image and generates prediction.	6	The system displays the result (e.g., "Benign" or "Malignant") with confidence % and store result in Database.
<b>Alternative Flow</b>			
3.1	Image is corrupted or unreadable.	3.2	The system displays: "Image format invalid or unreadable. Please upload a new image."

### 3.4.6 Add users:

**Table 7: Add users**

Name	Add Users		
Actors	Admin		
Summary	This use case describes the process of adding users to the SpotCancerAI platform.		
Pre-Conditions	The actor (Administrator) must have the necessary privileges and permissions to add users.  The actor must be logged in to their SpotCancerAI account.		
Post-Conditions	The new users are successfully added to the SpotCancerAI platform.		
Special Requirements	None		
Basic Flow			
Actor Action		System Response	
1	Actor logs in to their SpotCancerAI account with appropriate privileges (Administrator).		
2	Actor navigates to the user management or account		

	administration section of the SpotCancerAI platform.		
3	Actor selects the option to add users or create new accounts.	4	The system presents a form or interface to enter user details, such as email address, username, and other relevant information.
5	Actor fills in the required user information, ensuring the accuracy and completeness of the data.	6	The system validates the entered information, checking for any potential errors or conflicts (e.g., duplicate email addresses, invalid usernames).
7	Actor submits the form or clicks on the "Add User" button to initiate the user creation process.	8	The system processes the user creation request and generates a new user account.
		9	The system displays a confirmation message indicating that the user has been
			successfully added to the SpotCancerAI platform.
10	The new user receives an email or notification containing their login credentials and instructions for accessing the platform.		
<b>Alternative Flow</b>			
7.1	If there are any issues during the user creation process, such as validation errors or database constraints.	7.2	the system displays an error message and advises the actor to review and correct the provided information

### 3.4.7 Remove Users:

**Table 8:Remove Users**

<b>Name</b>	Remove user		
<b>Actors</b>	Admin		
<b>Summary</b>	This use case describes the process of removing users from the SpotCancerAI platform. The removal can only be performed by the administrator.		
<b>Pre-Conditions</b>	The administrator must have the necessary privileges and permissions to remove users.  The administrator must be logged in to their SpotCancerAI account.		
<b>Post-Conditions</b>	The selected users are successfully removed from the SpotCancerAI platform.		
<b>Special Requirements</b>	None		
<b>Basic Flow</b>			
<b>Actor Action</b>		<b>System Response</b>	
1	Administrator logs in to their SpotCancerAI account with appropriate privileges.		
2	Administrator navigates to the user management or account administration section of the SpotCancerAI platform.		
3	Administrator views the list of users on the platform.		
4	Administrator selects one or multiple users to be removed.		

5	Administrator confirms the selection and initiates the removal process.	6	The system prompts the administrator to confirm the removal action, ensuring they understand the consequences.
7	Administrator confirms the removal action.	8	The system processes the removal request and deletes the selected users' accounts from the SpotCancerAI platform.
		9	The system displays a confirmation message indicating that the users have been successfully removed.
<b>Alternative Flow</b>			
5.1	If there are any issues during the user removal process, such as database errors or system constraints.	5.2	The system displays an error message and advises the administrator to retry.

### 3.4.8 Check Results:

**Table 9:Check Results**

<b>Name</b>	Check Results
<b>Actors</b>	Admin, User, Dermatologists
<b>Summary</b>	The user can view the prediction results of previously analyzed skin images, including diagnosis, date, and confidence score.
<b>Pre-Conditions</b>	At least one image must have been analyzed by the system.
<b>Post-Conditions</b>	The prediction result is generated, stored, and displayed to the user.
<b>Special Requirements</b>	None
<b>Basic Flow</b>	

Actor Action		System Response	
1	The user logs into the system.	2	The system redirects to the dashboard/home page.
3	The user navigates to the "Results" or "History" section.	4	The system fetches the list of all analyzed images and their results.
5	The user selects a specific result to view details..	6	The system displays the result (diagnosis, confidence, date, and image preview).
<b>Alternative Flow</b>			
4.1	No results available for user.	4.2	The system displays: "No analysis results found. Please upload and analyze an image first."

### 3.4.9 View User History:

**Table 10: View User History**

<b>Name</b>	View User History
<b>Actors</b>	Admin, User, Dermatologists
<b>Summary</b>	The user can view a complete history of all uploaded images and their associated prediction results (if analyzed).
<b>Pre-Conditions</b>	The user must be logged into the system.
<b>Post-Conditions</b>	A list of all uploaded images and their statuses is displayed.
<b>Special Requirements</b>	None
<b>Basic Flow</b>	
<b>Actor Action</b>	<b>System Response</b>

1	The user clicks on “History” or “My Uploads.”	2	The system queries the database for all uploaded images associated with the user.
3	The system displays a list including image, upload date, status (Analyzed / Pending), and result (if available).		
4	The user selects a specific result to view details..	5	The system displays the result (diagnosis, confidence, date, and image preview).
<b>Alternative Flow</b>			
4.1	No results available for user.	4.2	The system displays: "No analysis results found.
5.1	Database fetch fails.	5.2	The system displays: "Unable to load history. Please try again later."

### 3.4.10 Give suggestions:

**Table 11: Give Suggestions**

<b>Name</b>	Give Suggestions
<b>Actors</b>	Dermatologists
<b>Summary</b>	After analyzing a skin image, the system or dermatologist provides suggestions such as recommended next steps, care tips, or referrals based on the prediction result.
<b>Pre-Conditions</b>	A prediction result must exist for the image.
<b>Post-Conditions</b>	Suggestions are displayed to the user and stored with the image record.

Special Requirements		None	
Basic Flow			
Actor Action		System Response	
1	The system completes image analysis and stores prediction.	2	System checks prediction result (e.g., Benign / Malignant)
3	System generates auto-suggestions (or dermatologist adds notes).	4	Suggestions are stored and linked with the corresponding result.
5	User views the result page.	6	System displays the diagnosis along with appropriate suggestions.
Alternative Flow			
4.1	No suggestion template available.	4.2	System shows a general message: “Consult a dermatologist for further guidance.”

### 3.4.11 Provide Feedbacks:

**Table 12:Provide Feedbacks**

<b>Name</b>	Provide feedback
<b>Actors</b>	User, Dermatologists
<b>Summary</b>	This use case describes the process of providing feedback on the SpotCancerAI.
<b>Pre-Conditions</b>	<input type="checkbox"/> The user must have a registered account on the SpotCancerAI platform. <input type="checkbox"/> The user must be logged in to their SpotCancerAI account.
<b>Post-Conditions</b>	The user's feedback is successfully submitted and received by the SpotCancerAI platform.



Special Requirements		None	
Basic Flow			
Actor Action		System Response	
1	User navigates to the SpotCancerAI website and logs in to their account.		
2	User locates the feedback section or finds the designated area for providing feedback.		
3	User clicks on the "Give Feedback" button.	4	The system presents a feedback form or text box.
5	User enters their feedback in the provided text area, providing specific details and information.		
6	User clicks on the "Submit" or "Send" button to submit the feedback.	7	The system captures and stores the user's feedback in the SpotCancerAI platform.
		8	The system displays a confirmation message indicating that the feedback has been successfully submitted.
Alternative Flow			


### 3.4.12 View FeedBacks:

**Table 13:FeedBacks**

Name	View feedback		
Actors	Admin		
Summary	The Admin can view all feedback submitted by users or Dermatologists, including comments on system performance, prediction quality, and feature suggestions.		
Pre-Conditions	The Admin must be logged in. At least one feedback must exist in the database.		
Post-Conditions	Feedbacks are retrieved and displayed with user details and timestamps.		
Special Requirements	None		
Basic Flow			
Actor Action		System Response	
1	Admin logs into the system.	2	System redirects to Admin Dashboard.
3	Admin clicks on the “View Feedbacks” module.	4	System fetches all feedback entries from the database.
5	Admin views feedback list, including: user name, role, date, feedback content, and associated prediction ID (if applicable).	6	System displays results with filter and sort options (e.g., by date, user type).
Alternative Flow			
4.1	System displays results with filter and sort options (e.g., by date, user type).	4.2	System displays: “No feedbacks have been submitted yet.”


## 3.5 GUI Graphical User Interfaces

### 3.5.1 Home Screen (Role Selection Interface)

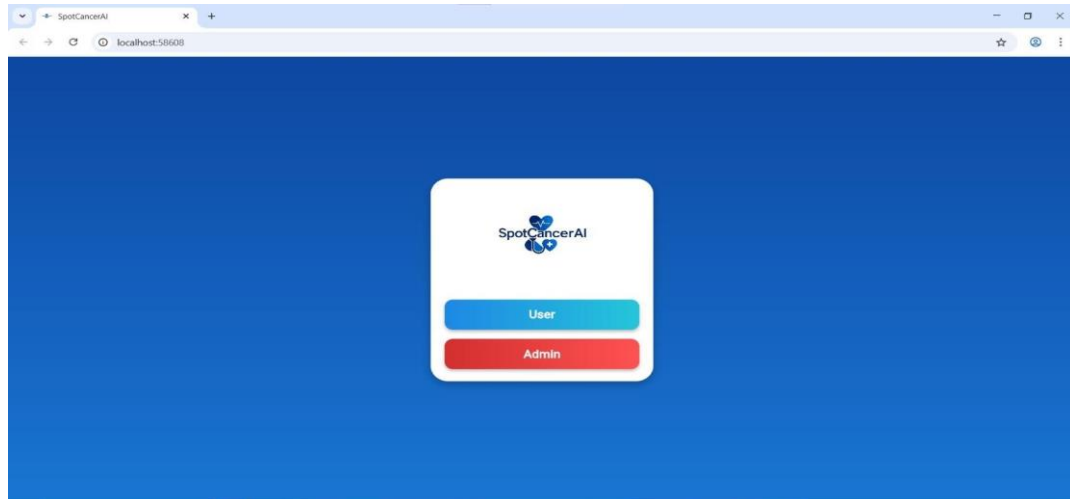


The Home Screen serves as the primary entry point of the SpotCancerAI system, presenting users with two main access options—User and Admin—to guide them toward the appropriate login module. The interface features a clean and modern design with a gradient blue background, while the SpotCancerAI logo is prominently displayed at the top of the centered card, creating a welcoming and professional appearance. This screen is designed to help users easily identify their role before proceeding, ensuring a clear separation between standard user functions and administrative access. The User button, styled with a blue gradient, directs individuals to the standard login page, whereas the Admin button, highlighted in a red gradient, leads to the administrator login panel, reinforcing the distinction between general and advanced system privileges.

### 3.5.2 User Login Screen

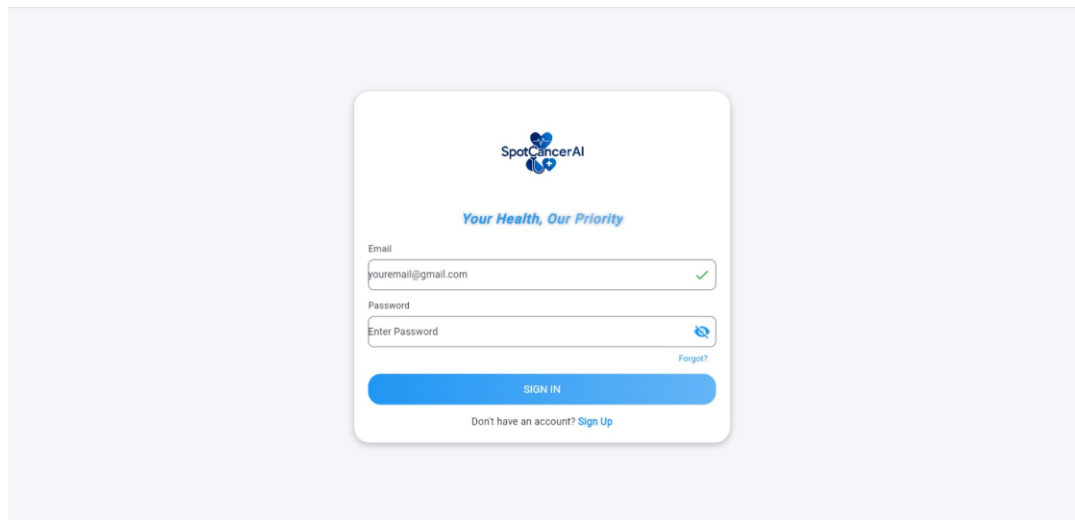

The User Login Screen provides a secure and user-friendly interface for patients or doctors to authenticate themselves before accessing the system's features. It incorporates a minimalistic layout with soft, light-toned colors that enhance readability and comfort for users. The screen includes fields for entering a registered email and password, with the password field offering a visibility toggle to improve usability. A clearly marked login button initiates the authentication process, while additional options such as "Forgot Password" allow users to recover access if necessary. The screen also includes a signup link, enabling new users to quickly create an account and begin using the system. Overall, the interface focuses on simplicity, accessibility, and a smooth login experience.


### 3.5.3 Admin Login Screen

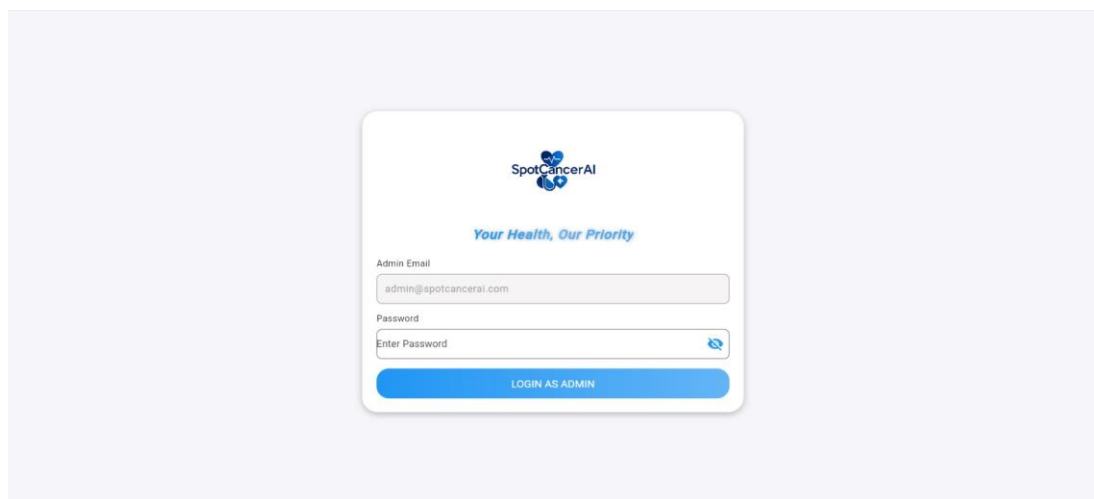

The Admin Login Screen is specifically designed for system administrators who manage backend operations, user permissions, and overall system security. It utilizes a bold red color theme that visually communicates the sensitive and restricted nature of administrative access. The screen includes dedicated fields for entering an admin email and password, both equipped with secure input formatting and a visibility toggle to reduce login errors.

The “Login as Admin” button validates credentials and grants access to the administrative dashboard. Additionally, a distinct security notice is displayed to remind users that the section is intended only for authorized personnel. The overall design emphasizes security, clarity, and controlled access to ensure proper handling of system-level functionalities.


### 3.5.4 Dashboard Screen (Main User Interface)

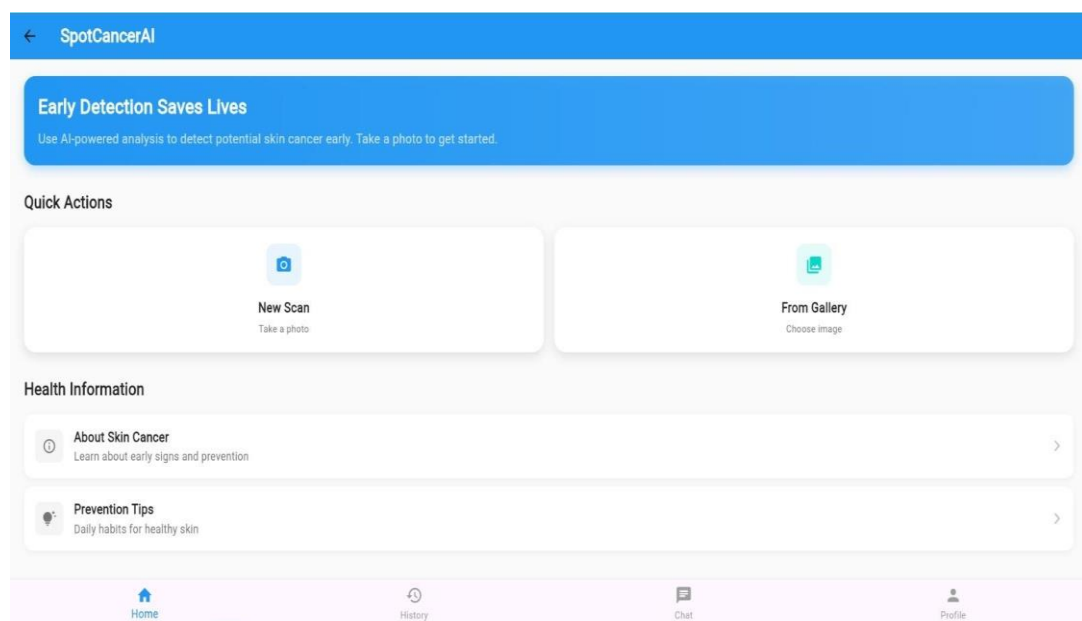

The Dashboard Screen serves as the central hub of the SpotCancerAI application, providing users with seamless access to all major functionalities in a clean and intuitive layout. From this screen, users can upload a dermoscopic image directly from their device or capture a new image using the built-in camera feature for instant analysis. The dashboard also includes a dedicated section for health information, where users can learn about different skin lesion types, prevention

methods, and early detection guidelines. A "Tips" panel offers expert-recommended practices for maintaining healthy skin and identifying warning signs of melanoma. Additionally, users can view their complete analysis history, allowing them to track previous scans and monitor any changes over time. The screen also integrates a contact option for reaching out to support or medical professionals when assistance is needed. A built-in chat feature enables real-time communication with healthcare experts or system support. Finally, a profile section allows users to manage their personal information, account settings, and preferences. Overall, the dashboard is designed for easy navigation, delivering a smooth and informative user experience tailored to both patients and medical practitioners.

### 3.5.5 Analysis Result Screen

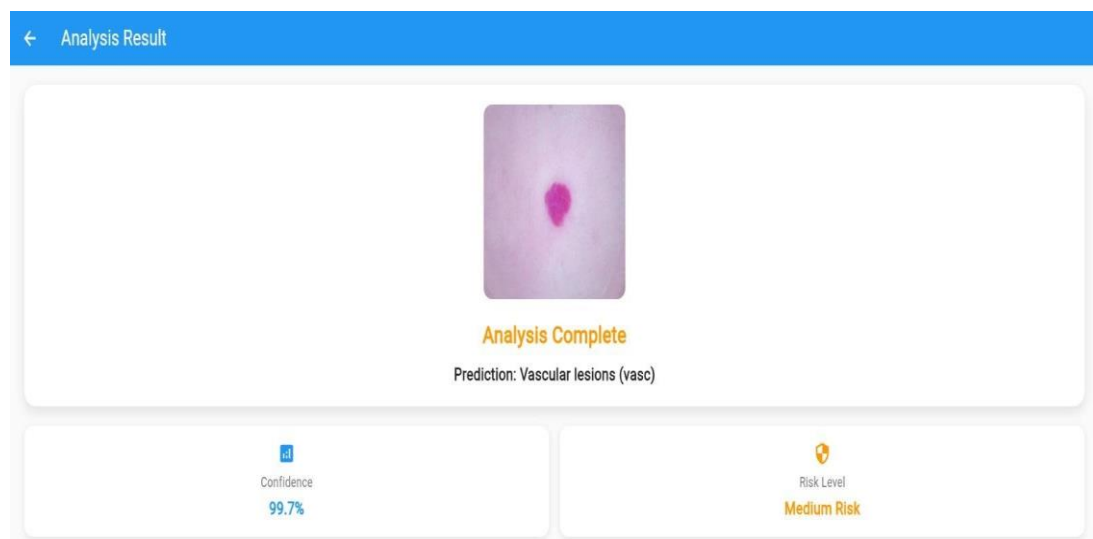

The Analysis Result Screen presents a clear and comprehensive summary of the outcome generated by the SpotCancerAI model after processing the uploaded or captured skin image. This interface displays the predicted class of the skin lesion—such as melanoma, nevus, or keratosis—allowing users and clinicians to immediately understand the model's diagnostic interpretation. Along with the predicted category, the screen provides the confidence score, which quantifies how certain the model is about its prediction, offering additional insight into decision reliability. The system also calculates and

displays a risk level (low, moderate, or high), helping users gauge the urgency of seeking medical attention. The results are presented in a visually structured format, often accompanied by color-coded indicators to enhance readability and support quick evaluation. Overall, the Analysis Result Screen is designed to deliver accurate, interpretable, and user-friendly diagnostic feedback that supports informed medical decision-making.

## **Chapter 4: Implementation and Test Cases**



## **4.1. Introduction**

The implementation of the proposed skin cancer classification system involves the development of a deep learning-based pipeline designed to accurately identify and categorize various types of skin cancer lesions. The system processes dermatological images through a sequence of well-defined stages, beginning with the detection of skin regions and culminating in the classification of skin cancer into seven distinct categories. These categories include Actinic Keratoses (akiec), Basal Cell Carcinoma (bcc), Benign Keratosis-like Lesions (bkl), Dermatofibroma (df), Melanoma (mel), Nevus (nv), and Vascular Lesions (vasc).

This section presents a comprehensive overview of the proposed methodology, encompassing dataset preparation, pipeline architecture, and model training and evaluation.

## **4.2. Implementation**

### **4.2.1. Proposed Framework**

We propose a deep learning-based skin cancer classification pipeline using the HAM10000 dataset, which includes seven categories of skin lesions. The system employs an EfficientNetB5-based Convolutional Neural Network (CNN) to perform precise multi-class classification across Actinic Keratoses (akiec), Basal Cell Carcinoma (bcc), Benign Keratosis-like Lesions (bkl), Dermatofibroma (df), Melanoma (mel), Nevus (nv), and Vascular Lesions (vasc). The model leverages transfer learning and extensive data preprocessing to enhance feature extraction and classification accuracy. This approach provides a reliable, efficient, and scalable framework for automated skin cancer diagnosis.

#### **4.2.1.1. Pipeline Overview**

##### **4.2.1.1.1. Multi-Class Classification with EfficientNet**

The proposed system employs a seven-class skin cancer classification model based on the EfficientNet architecture, trained on the HAM10000 dataset. Since the dataset exclusively contains dermoscopic images of skin lesions, no preliminary filtering is required. The model performs multi-class classification across seven categories: Actinic Keratoses (*akiec*), Basal Cell Carcinoma (*bcc*), Benign Keratosis-like Lesions (*bkl*), Dermatofibroma (*df*), Melanoma (*mel*), Nevus (*nv*), and Vascular Lesions (*vasc*).

To address the dataset's class imbalance, techniques such as data augmentation and class weighting are applied during training. The EfficientNet model is fine-tuned using transfer learning to leverage pre-trained ImageNet features for improved accuracy and generalization. This stage serves as the core of the pipeline, enabling precise and efficient classification of skin cancer lesions based on dermoscopic imagery.

#### 4.2.1.1.2. Multi-Class Classification with Ensemble Model

The second stage of the proposed pipeline employs a multi-class skin cancer classification model trained on the HAM10000 dataset, which includes seven diagnostic categories: *Actinic Keratoses (akiec)*, *Basal Cell Carcinoma (bcc)*, *Benign Keratosis-like Lesions (bkl)*, *Dermatofibroma (df)*, *Melanoma (mel)*, *Nevus (nv)*, and *Vascular Lesions (vasc)*. To enhance performance, an ensemble model combining EfficientNet-B5 and EfficientNet-B1 architectures was implemented. The ensemble approach aimed to leverage the complementary strengths of both networks to improve generalization and classification accuracy.

Although class balancing techniques such as data augmentation and class weighting were applied during training, the ensemble model did not yield significant improvements compared to the individual models. This performance limitation may be attributed to overlapping feature representations between the two EfficientNet variants or insufficient model diversity. Future enhancements could involve combining more heterogeneous architectures, employing

Advanced ensemble fusion strategies (e.g., stacking or weighted averaging), and refining data balancing methods to better address class distribution disparities within the HAM10000 dataset.

#### **4.2.1.1.3 Performance of Inception and ResNet Models on HAM10000**

In the third stage, Inception and ResNet architectures were individually trained on the HAM10000 dataset for seven-class skin cancer classification. However, both models demonstrated poor performance and failed to achieve satisfactory accuracy or generalization. This limited performance can be attributed to the high class imbalance inherent in the dataset, as well as the visual similarity among certain lesion types, which makes fine-grained classification challenging.

Moreover, both architectures struggled to extract discriminative features due to variations in image quality, lighting, and lesion texture. Despite applying standard preprocessing, data augmentation, and class balancing techniques, the models exhibited overfitting and low sensitivity for minority classes such as Actinic Keratoses (akiec) and Dermatofibroma (df). The results suggest that Inception and ResNet are not well-suited for this specific dataset without significant architectural adaptation or specialized preprocessing.

#### **4.2.1.1.4 Final Multi-Class Skin Cancer Classification Using**

##### **EfficientNetB5**

In the final stage, the EfficientNetB5 architecture was employed for seven-class skin cancer classification using the HAM10000 dataset. To address the challenges of class imbalance and improve model robustness, data augmentation and focal loss were applied during training. Random augmentation techniques, including rotation, flipping, zooming, and brightness adjustments, were utilized to increase dataset diversity and reduce overfitting.

The model demonstrated strong performance, achieving over 90% training accuracy and 84% validation accuracy, indicating effective generalization across unseen data. The integration of focal loss helped the model focus more on minority classes, enhancing sensitivity and reducing bias toward dominant categories such as Nevus (nv).

Overall, this stage produced the best results among all configurations, confirming that the combination of EfficientNetB5, focal loss, and random augmentation forms a highly effective framework for accurate and balanced skin cancer classification using the HAM10000 dataset.

The pipeline's architecture is illustrated in the Figure, which depicts the flow from input image processing through skin detection to disease classification. This visual representation highlights the sequential nature of the framework and the role of each model in achieving the overall objective of accurate skin disease diagnosis.


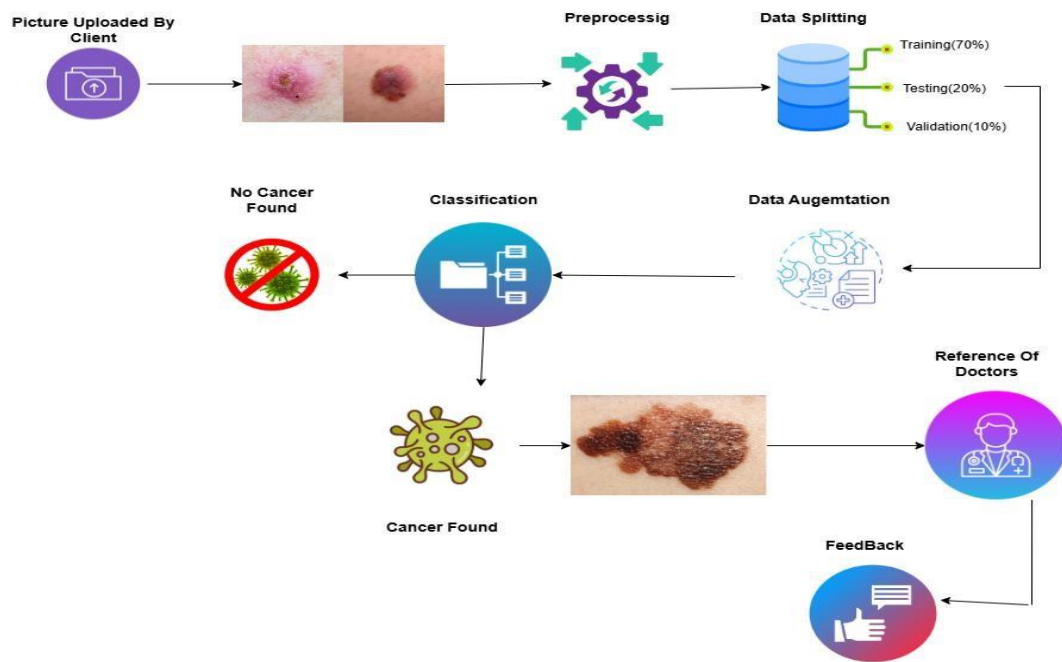


Figure 3:Skin Cancer Classification

#### 4.2.2. Dataset

We utilized the HAM10000 dataset for our CNN-based skin disease classification pipeline. This dataset contains 10,015 dermatoscopic images across seven disease categories, enabling both robust preprocessing and accurate multi-class classification. The distribution spanning melanocytic nevi (6,705), melanoma (1,113), benign keratosis (1,099), basal cell carcinoma (514), actinic keratoses (327), vascular lesions (142), and dermatofibroma (115) informed our class balancing strategy to ensure reliable model training and evaluation.

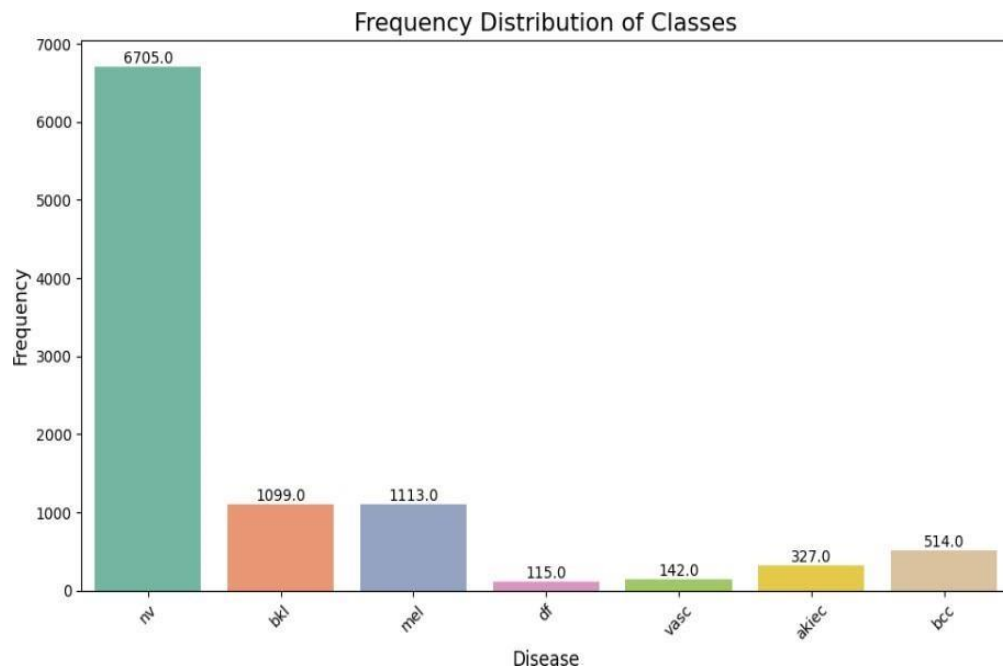


**Figure 4:Distribution of Classes**

#### **4.2.2.1 Solving Data Imbalance**

The HAM10000 dataset suffers from severe class imbalance, with melanocytic nevi dominating (67%) while rare classes like dermatofibroma are underrepresented (1.1%). This imbalance would bias the model toward majority classes.


We solved this through data augmentation applying rotations, flips, and brightness adjustments—with more aggressive augmentation for minority classes. Combined with class- weighted loss functions that assigned higher penalties for misclassifying rare diseases, this approach ensured balanced learning across all categories and improved detection of both common and rare conditions.

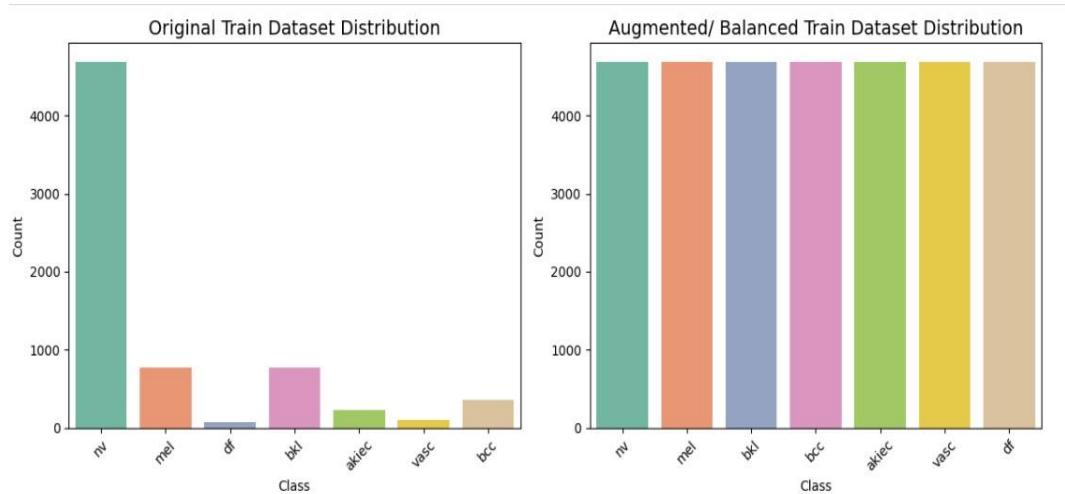

**Figure 5:Balancing of Classes**

### 4.2.3 Preprocessing Image

The preprocessing pipeline standardizes dermatological images through systematic hair removal using morphological operations and exemplar-based inpainting, followed by dimensional normalization to  $456 \times 456$  pixels using Lanczos interpolation. Color consistency is maintained through RGB conversion, while adaptive thresholding ensures precise detection of occluding structures. This comprehensive preprocessing framework preserves critical diagnostic features and ensures input homogeneity for subsequent convolutional neural network analysis, significantly enhancing model reliability and classification accuracy.

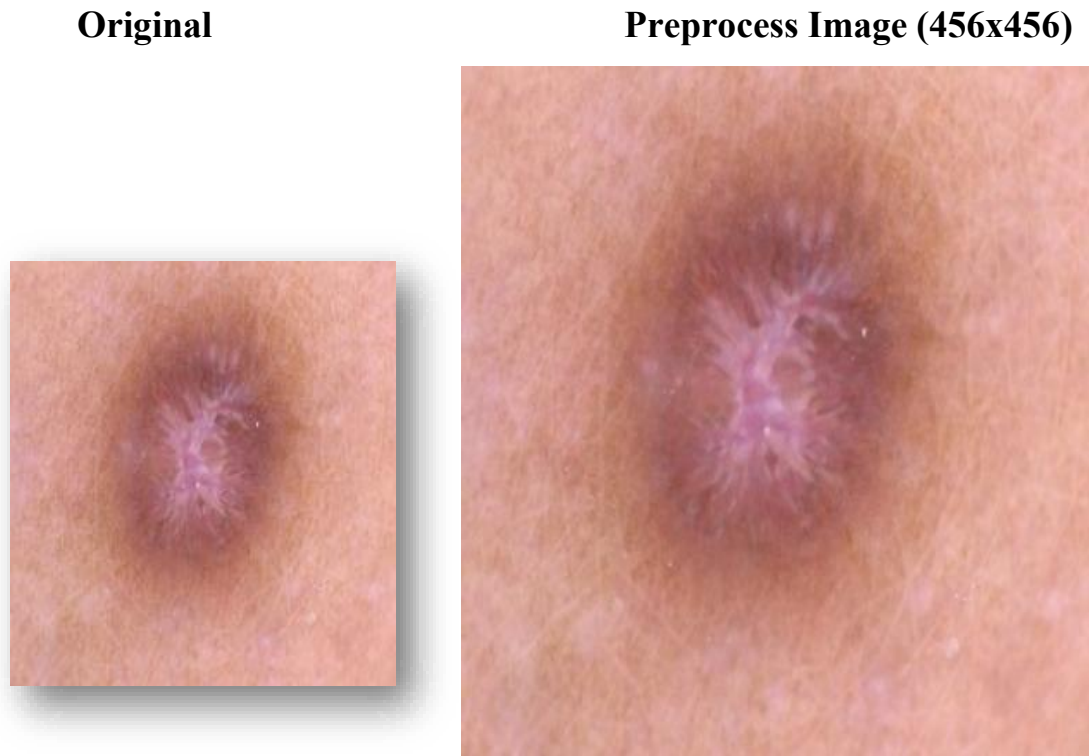

**Figure 6:Image Preprocessing**

#### 4.2.4 Dataset Splitting

The HAM10000 dataset was partitioned into standardized training, validation, and test subsets to facilitate robust model development and evaluation. The dataset was divided into 7,010 images (70.00%) for training, 1,001 images (10.00%) for validation, and 2,004 images (20.01%) for testing.

##### **Before Balance:**

**Table 14: Dataset Before Balancing**

Split Type	Number of Images	Percentage
<b>Train</b>	7,010	70.00%
<b>Validation</b>	1,001	10.00%
<b>Test</b>	2,004	20.00%
<b>Total</b>	<b>10,015</b>	<b>100.00%</b>

## After Balance:

**Table 15: Dataset After Balancing**

Split Type	Number of Images	Percentage	Notes
<b>Train</b>	32,851	91.60%	After augmentation & balancing
<b>Validation</b>	1,001	2.79%	Original samples only
<b>Test</b>	2,004	5.59%	Original samples only
<b>Total</b>	<b>35,856</b>	<b>100%</b>	After balancing procedures

### 4.2.5 CNN Architecture

The convolutional neural network architecture developed for dermatological image classification employs a sophisticated transfer learning framework centered on EfficientNetB5 as the core feature extractor. This backbone network, pre-trained on the extensive ImageNet dataset, provides a robust foundation for visual pattern recognition, which is particularly advantageous for medical imaging tasks with limited annotated data. The architectural implementation begins with strategic fine-tuning methodology, where the initial 200 layers remain frozen to preserve generic feature detection capabilities, while subsequent layers are made trainable to adapt the network to domain-specific dermatological characteristics and lesion patterns.

The custom classification head demonstrates a carefully engineered progressive dimensionality reduction approach, transitioning through 512, 256, and 128 fully-connected units. Each dense layer incorporates Swish activation functions, recognized for superior performance in deep networks compared to traditional ReLU, accompanied by L2 regularization ( $\lambda=1e-4$ ) to constrain weight magnitudes and prevent overfitting. Batch normalization layers are systematically integrated after each dense layer to stabilize training dynamics and accelerate convergence through internal covariance shift reduction. A descending dropout strategy is implemented with rates of 0.5, 0.3, and 0.2 respectively, providing progressively



stronger regularization in earlier layers while allowing more precise feature learning in deeper layers.

The optimization framework addresses the critical challenge of class imbalance through a customized categorical focal loss function with hyperparameters  $\gamma=3.0$  and  $\alpha=0.5$ , which dynamically scales the loss to focus learning on misclassified examples and balance class contributions. Model compilation utilizes the Adam optimizer with a conservative learning rate of  $1e-5$  to ensure stable gradient updates during fine-tuning. Comprehensive evaluation metrics include standard categorical accuracy, area under the ROC curve (AUC) for diagnostic capability assessment, and top-3 categorical accuracy to capture clinical relevance where multiple differential diagnoses may be considered. This multi-faceted architectural design ensures robust feature extraction, effective regularization, and clinically meaningful performance evaluation for automated skin disease classification.

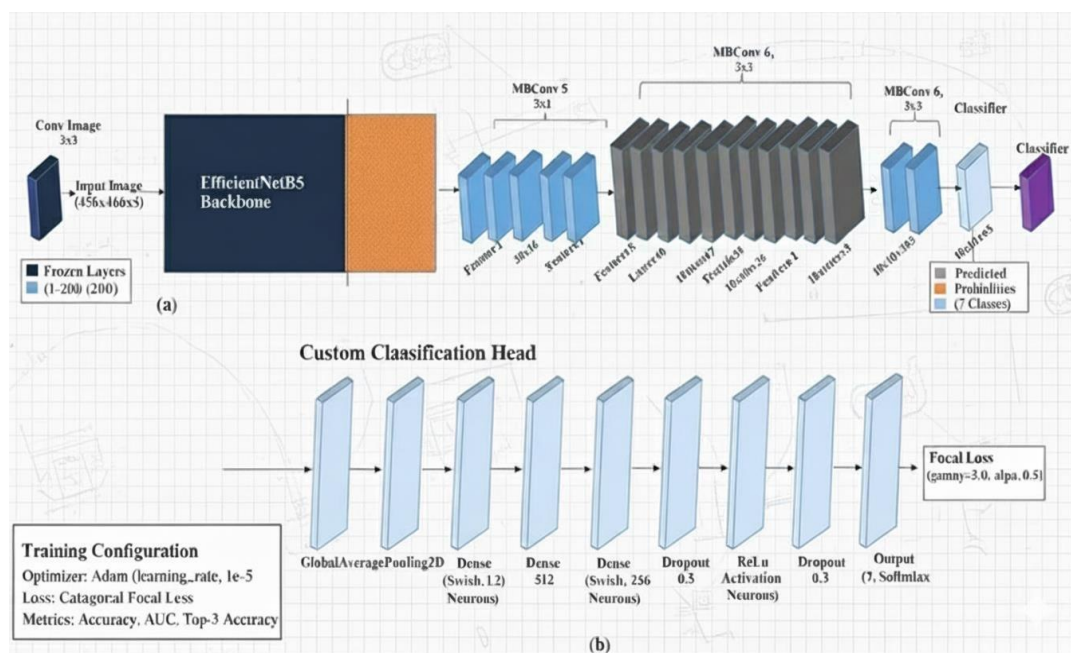


Figure 7: CNN Architecture

#### 4.2.6 Activation and Loss Function Activation Function (Softmax):

Used in the output layer to convert logits into class probabilities.

$$\sigma(z_i) = \frac{e^{z_i}}{\sum_{j=1}^K e^{z_j}}$$

It ensures all output values are between 0 and 1 and sum to 1, suitable for multi-class classification. **Loss Function (Categorical Focal Loss):**

Focuses training on hard, misclassified examples by reducing the weight of well-classified ones.

$$FL(p_t) = -\alpha(1 - p_t)^\gamma \log(p_t)$$

Here,  $p_t$  is the predicted probability for the true class,  $\alpha$  balances class weights, and  $\gamma$  controls the focusing strength.

#### 4.2.7 Model Training

The model training was performed on **Kaggle's GPU environment**, which provided the necessary computational power for deep learning operations. The proposed architecture utilized **EfficientNetB5** as the base model with ImageNet pre-trained weights for transfer learning, enabling efficient feature extraction from skin lesion images. Custom fully connected layers with **Swish activations**, **Batch Normalization**, and **Dropout regularization** were added to enhance learning stability and prevent overfitting. The model was fine-tuned by unfreezing the last layers of EfficientNetB5 to adapt to the HAM10000 dataset. The training process employed the **Adam optimizer** with a **learning rate of 1e-5**, chosen for its adaptive gradient handling and stable convergence. The **categorical focal loss** function was used to focus learning on harder, misclassified samples, improving robustness against class imbalance. Performance metrics included **accuracy**, **AUC**, and **Top-3 categorical accuracy**, allowing comprehensive evaluation of the model's classification capability across seven skin disease categories.

#### 4.2.8 Model Evaluation

The evaluation of the trained model was conducted using standard performance metrics, including **Precision**, **Recall**, **F1-Score**, and the **Confusion Matrix**, comprehensively assess its classification effectiveness. **Precision** measures the proportion of correctly predicted positive samples among all predicted positives, reflecting the model’s accuracy in identifying true cases. **Recall** (or Sensitivity) evaluates the model’s ability to correctly detect all actual positive cases, ensuring minimal false negatives. The **F1-Score**, being the harmonic mean of Precision and Recall, provides a balanced assessment of the model’s accuracy and robustness, particularly in handling class imbalance. Additionally, the **Confusion Matrix** visually represented the distribution of true and false classifications across all seven skin disease categories, enabling detailed insight into model performance and misclassification trends. Together, these metrics confirmed the model’s reliability and strong discriminative capability in accurately identifying various types of skin lesions.

4.2.9 Conclusion

The implementation of the proposed pipeline demonstrates a systematic approach to handling complex datasets and achieving effective classification. Despite challenges such as dataset imbalance and noise, the system effectively identifies and classifies skin diseases in a two-stage process. Continuous dataset improvements and model optimizations are expected to enhance performance further.

4.3 Test Case\_1

Table 16: Test Case1

SpotCaancerAI/Patient/Signup Module			
FYP II Documentation Section 4.3.1			
Test Case ID:	TA-01	Test Date:	3-3-2025
Test case Version:	V1.0	Use Case Reference(s):	UC-Patient Signup
Revision History:	NILL		
Objective	Testing the Signup module for the Patients.		
Product/Ver/Module:	SpotCancerAI/Patient/Signup		
Environment:	PC/ Browser/internet connectivity		

<b>Assumptions:</b>		<i>Patient Signup for the first time on system.</i>
<b>Pre-Requisite:</b>		<i>Have access the SpotCancerAI Signup page of the patient.</i>
<b>Step No.</b>	<b>Execution description</b>	<b>Procedure result</b>
<b>01</b>	<i>Enter detail in the relevant fields with correct formatting.</i> <i>(Name: {A.... Z, a.... z}.</i> <i>Email: Contains '@' along with</i> <i>{{a, b, c.... z}, {A, B, C,...Z},</i> <i>{1,2,3,...}}}</i>	<i>System registers the user with pop-up indicating Signup successfully.</i>
	<i>before '@' {.,_-} before '@'</i> <i>Ending with '.' And after those</i> <i>some alphabetic letters</i> <i>Password length &gt;= 6</i> <i>Contact: {1,2,3....}).</i>	
<b>02</b>	<i>Enter invalid format of field</i> <i>"Name".</i> <i>{1,2,3....} {#, @ and other</i> <i>special characters} OR</i> <i>"Contact" Multiple or Special</i> <i>characters.</i> <i>OR the field "Password" Length</i> <i>&lt; 6. OR the field "Email" Not</i> <i>contains '@' Multiple '@'</i> <i>characters Not</i> <i>ending at '.' Without letters.</i>	<i>System indicates error message respectively "incorrect name format", "incorrect contact", "incorrect password", "incorrect email". And doesn't register the user.</i>
<b>04</b>	<i>Enter duplicate email which</i> <i>already</i> <i>registered in system.</i>	<i>System indicate error "email already exist"</i>
<b>Comments:</b> <b>Only Valid Formatting of the fields are accepted.</b> <b>Duplicate emails not accepted.</b>		
<i>Passed</i> <i>Failed</i> <i>Not Executed</i>		

### 4.3.1 Test Case\_2

Table 17: Test Case 2

SpotCancerAI/Admin/Login Module			
FYP II Documentation Section 4.3.2			
<b>Test Case ID:</b>	TA-02	<b>Test Date:</b>	3-7-2025
<b>Test case Version:</b>	V1.0	<b>Use Case Reference(s):</b>	UC-Admin Login
<b>Revision History:</b>	NILL		
<b>Objective</b>	Testing the login module for the Admin.		
<b>Product/Ver/Module:</b>	SpotCancerAI/Admin/Login		
<b>Environment:</b>	PC/ Browser		
<b>Assumptions:</b>	Admin logs into the system using a specific email and password.		
<b>Pre-Requisite:</b>	Have access the SpotCancerAI Admin Login Page.		
Step No.	Execution description	Procedure result	
01	<p>Enter valid login credentials in the respective fields with correct formatting.</p> <p>(Email: Contains '@' along with {a, b, c.... z}, {A, B, C,....Z}, {1,2,3,...})</p> <p>Before '@': {.,_-} allowed.</p> <p>Ending with '.' and followed by alphabetic letters.</p> <p>Password: Must match the specific password assigned to the admin account and have a length <math>\geq 6</math>.)</p>	<p>System verifies the credentials and successfully logs in the admin, displaying the <b>Admin Dashboard</b>.</p>	

02	<p><i>Enter invalid or incorrect credentials in either the “Email” or “Password” field.</i></p> <p><i>Examples:</i></p> <ul style="list-style-type: none"> <li>• <i>Email not containing ‘@’ or containing multiple ‘@’ characters.</i></li> <li>• <i>Email not ending with a valid domain (e.g., “.com”).</i></li> <li>• <i>Password incorrect or length &lt; 6.</i></li> </ul>	<p>System displays an appropriate error message: <b>“Incorrect email or password.”</b> The system does <b>not</b> allow access to the admin dashboard.</p>
<p align="center"><b>Comments:</b></p> <p align="center"><b>Only Valid Formatting of the fields are accepted.</b></p>		
	<div>Passed</div> <div></div>	<div>Failed</div> <div></div> <div>Not Executed</div> <div></div>

### 4.3.2 Test Case\_3

Table 18:Test Case 3

SpotCancerAI/Patient/Login Module			
FYP II Documentation Section 4.3.3			
Test Case ID:	TA-03	Test Date:	3-10-2025
Test case Version:	V1.0	Use Case Reference(s):	UC-Patient Login

<b>Revision History:</b>		<i>NILL</i>
<b>Objective</b>		<i>Testing the Login module for the Patients.</i>
<b>Product/Ver/Module:</b>		<i>SpotCancerAI/Patient/Login</i>
<b>Environment:</b>		<i>PC/ Browser</i>
<b>Assumptions:</b>		<i>User data exists in the system database.</i>
<b>Pre-Requisite:</b>		<i>Must Signup/register in the system before.</i>
<b>Step No.</b>	<b>Execution description</b>	<b>Procedure result</b>
<b>01</b>	<i>Enter correct Email and password.</i>	<i>System matches the credentials from the database and indicates the pop-up message that “login successfully”</i>
<b>02</b>	<i>Enter invalid email and correct password.</i>	<i>System matches the credentials from database indicates error “incorrect email”</i>
<b>03</b>	<i>Enter valid email and incorrect password.</i>	<i>System indicates error “incorrect password”</i>
<b>04</b>	<i>Enter invalid email and invalid password.</i>	<i>System indicates error “incorrect email or password”.</i>
<b>05</b>	<i>Enter credentials not registered in the system. OR enter name and password field remains empty.</i>	<i>System indicates error “user doesn’t exist”</i>
<b>Comments:</b>		
<b>Only valid user can login to the system who registered himself prior.</b>		
<div> <div>Passed</div> <div>Failed</div> <div>Not Executed</div> </div>		

### 4.3.3 Test Case\_4

Table 19: Test Case 4

<b>SpotCancerAI/Patient/Login Module</b>			
<b>FYP II Documentation Section 4.3.4</b>			
<b>Test Case ID:</b>	<i>TA-05</i>	<b>Test Date:</b>	<i>3-24-2025</i>

<b>Test case Version:</b>	<i>V1.0</i>	<b>Use Case Reference(s):</b>	<i>UC-Patient Change Password</i>
<b>Revision History:</b>	<i>NILL</i>		
<b>Objective</b>	<i>Testing the Change Password module for the Patient.</i>		
<b>Product/Ver/Module:</b>	<i>SpotCancerAI/Patient/Change Password</i>		
<b>Environment:</b>	<i>PC/ Browser</i>		
<b>Assumptions:</b>	<i>User is logged in and has a registered email in the system database.</i>		
<b>Pre-Requisite:</b>	<i>Patient must be registered and logged into the system.</i>		
<b>Step No.</b>	<b>Execution description</b>	<b>Procedure result</b>	
<b>01</b>	Enter correct current password, valid new password, and matching confirm new password.	System validates credentials and displays a pop-up message “Password changed successfully.”	
<b>02</b>	Enter incorrect current password, valid new password, and matching confirm new password.	System indicates error “Incorrect current password.”	
<b>03</b>	Enter correct current password, new password that does not meet requirements (e.g., too short), and matching confirm new password.	System indicates error “New password does not meet requirements.”	
<b>04</b>	Enter correct current password, valid new password, and non-matching confirm new password.	System indicates error “Confirm password does not match.”	
<b>05</b>	Leave any field (current password, new password, or confirm new password) empty.	System indicates error “All fields are required.”	



<b>Comments:</b>
Only valid user can login to the system who registered himself prior.
<b>Passed</b> <b>Failed</b> <b>Not Executed</b>

### 4.3.4 Test Case\_5

Table 20: Test Case 5

SpotCancerAI/Patient/Login Module			
FYP II Documentation Section 4.3.5			
<b>Test Case ID:</b>	TA-05	<b>Test Date:</b>	4-1-2025
<b>Test case Version:</b>	V1.0	<b>Use Case Reference(s):</b>	UC-Admin Account Management
<b>Revision History:</b>	NILL		
<b>Objective</b>	Testing the Account Management in the Admin Dashboard for activating, deactivating, and deleting patient accounts.		
<b>Product/Ver/Module:</b>	SpotCancerAI/Admin/Dashboard		
<b>Environment:</b>	PC/ Browser		
<b>Assumptions:</b>	Admin is logged in with appropriate permissions to manage accounts; patient accounts exist in the system.		
<b>Pre-Requisite:</b>	Admin logged in, and have access to the Dashboard.		
Step No.	Execution description	Procedure result	
01	Log in as an admin, navigate to the Dashboard, select an inactive patient account, and activate it.	System activates the patient account; the patient can now access their account.	
02	Log in as an admin, navigate to the Dashboard, select an active patient account, and deactivate it; attempt to log in as the patient.	System deactivates the patient account; the patient cannot access the dashboard (login fails with error).	

<b>03</b>	Log in as an admin, navigate to the Dashboard, select an active patient account, and delete it; attempt to log in as the patient.	System deletes the patient account; the patient cannot access the dashboard (login fails with error).
<b>Comments:</b>  <b>Admin can activate, deactivate, and delete patient accounts. Deactivated or deleted accounts prevent users from accessing their dashboards. Only existing patient accounts can be managed.</b>		
<i>Passed</i> <i>Failed</i> <i>Not Executed</i>		

### 4.3.5 Summary

In this section, we summarize the common attributes and metrics used to evaluate the test cases conducted for the **SpotCancerAI** platform, ensuring the reliability and functionality of its modules, such as **Patient Signup, Patient Login, Admin Login, Account Management, and Change Password**. The test metrics provide a standardized framework to measure the success, efficiency, and coverage of the testing process. These metrics are derived from the attributes shared across all test cases and are essential for assessing the system's performance in real-world scenarios.

The following are the key test metrics and their common attributes:

1. **Test Case Objective:** Each test case is designed with a clear objective to validate specific functionalities of the system. For instance, the objective may include testing patient registration, login verification, account management operations, or password updates, ensuring alignment with functional requirements outlined in the use case references.
2. **Test Execution Status:** The outcome of each test case is categorized as "Passed," "Failed," or "Not Executed." This metric indicates whether the system behaves as expected under the defined conditions. For example, in account management test cases, all steps passed, confirming correct functionality of patient account operations.
3. **Test Coverage:** This metric evaluates the extent to which the system's functionalities are tested. Each test case targets a specific module (e.g., Patient Signup, Admin Login) and verifies critical operations, ensuring

comprehensive coverage of the system’s features.

- 4. **Pre-requisites and Assumptions:** Common attributes include pre-conditions required for executing the test, such as user authentication (e.g., admin logged in with appropriate permissions) and system state (e.g., existing patient accounts). These ensure the test environment is consistent and replicable.
- 5. **Procedure and Result Validation:** Each test case includes a step-by-step execution procedure and corresponding results. The results are validated against expected outcomes, such as successful patient registration or failure to log in after deactivation. This metric ensures that the system’s behavior aligns with the intended design.
- 6. **Error Handling and Comments:** This attribute captures any issues, observations, or system limitations encountered during testing. Comments provide insights into the system’s performance, such as restrictions on managing non-existing accounts or error messages for failed logins.
- 7. **Test Environment:** The test environment, including platform (e.g., PC/Browser) and module (e.g., SpotCancerAI/Admin/Dashboard), is consistently documented to ensure tests are conducted under standardized conditions. This metric supports reproducibility and scalability of the testing process.

These metrics collectively provide a quantitative and qualitative assessment of the system’s functionality, enabling stakeholders to gauge the reliability and robustness of the **SpotCancerAI** platform. By analyzing the outcomes of test cases across modules such as **Patient Signup, Patient Login, Admin Login, Account Management,** and **Change Password**, we ensure that the system meets user requirements and performs effectively in managing critical operations.

4.4 Test Case Metric

Table 21: Test Case Metric

Metric:	Purpose
---------	---------

<b>Number of Test Cases:</b>	5
<b>Number of Test Cases Passed:</b>	5
<b>Number of Test Cases Failed:</b>	0
<b>Test Case Defect Density:</b>	$(0 \times 100) / 5 = 0\%$
<b>Test Case Effectiveness:</b>	If all defects were found via test cases: $(5/5 \times 100) = 100\%$
<b>Traceability Matrix:</b>	All 5 test cases map directly to defined use cases for modules Patient Signup, Patient Login, Admin Login, Account Management, and Change Password. Traceability maintained.

## **Chapter 5: Experimental Results and Analysis**

## 5.1 Introduction


This chapter presents a comprehensive evaluation of the experimental results and analysis for the skin disease classification system using the HAM10000 dataset. We detail the experiments conducted with our proposed EfficientNetB5-based architecture, specifically designed to address the challenges of dermatological image classification, including class imbalance and fine-grained feature extraction. The model's performance is systematically compared against several state-of-the-art convolutional neural network architectures to demonstrate its efficacy and competitive advantages. Through rigorous analysis of multiple performance metrics and visualization techniques, we provide insights into the model's strengths, limitations, and potential areas for future improvement in automated skin disease diagnosis.

## 5.2 Experiments

We conducted comprehensive experiments on the HAM10000 dataset to evaluate our proposed EfficientNetB5-based model for multi-class skin disease classification. Our architecture incorporates strategic fine-tuning, focal loss for class imbalance, and advanced regularization techniques. For benchmarking, we compared against three state-of-the-art models: **InceptionV3** for multi-scale feature extraction, **ResNet50** for deep residual learning, and standard **EfficientNet** variants for optimized scaling. All models were evaluated under identical conditions using the same preprocessing pipeline and balancing strategies to ensure fair comparison. This section presents detailed performance metrics and comparative analysis.

### 5.2.1 EfficientNetB5 Model Performance (OverSample)

The performance of the EfficientNetB5-based skin lesion classification model demonstrates strong learning and generalization capabilities. From the training curves, the left plot shows the loss convergence for both training and validation sets over 21 epochs. The training loss decreases consistently from a high initial value, reaching approximately 0.16 by the final epoch, indicating that the model is effectively learning patterns from the training data. The validation loss decreases initially but stabilizes around 0.22, which suggests that the model generalizes reasonably well without significant overfitting. On the right, the accuracy curves show that the training accuracy steadily increases, reaching around 92%, while the validation accuracy stabilizes at approximately 84%. This gap indicates some degree of overfitting, but overall, the model maintains strong predictive performance on unseen data.

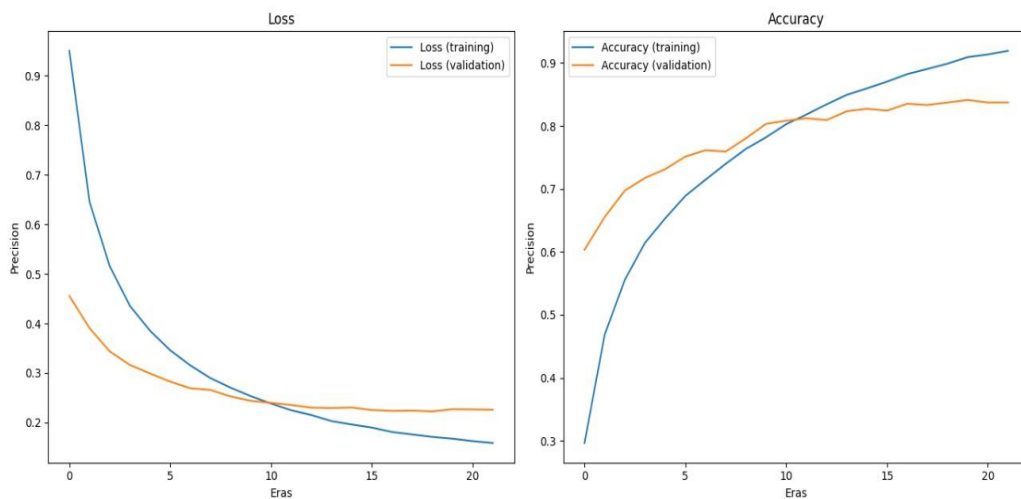

**Figure 8: Efficient Net B5 Model Graph**

## Classification Report:

The classification report further confirms the model's performance across different skin lesion classes. The model achieves the highest performance on the nv class (normal skin) with a precision of 0.93, recall of 0.91, and F1-score of 0.92, reflecting its ability to correctly identify the majority class. Moderate performance is observed for other classes such as bcc (precision 0.90, recall 0.76) and vasc (precision 0.89, recall 0.83). Lower performance is noted for challenging and underrepresented classes like mel

(precision 0.55, recall 0.64) and akiec (precision 0.73, recall 0.69), which indicates room for improvement, potentially through additional data augmentation or class-balancing techniques. Overall, the weighted F1-score of

0.84 and accuracy of 84% demonstrate that the model performs robustly across all classes, making it suitable for practical skin lesion classification while highlighting areas where minority classes could benefit from further enhancement.

**Table 22:Classification Report**

<b>Class</b>	<b>Precision</b>	<b>Recall</b>	<b>F1-Score</b>	<b>Support</b>
<b>akiec</b>	0.73	0.69	0.71	65
<b>bcc</b>	0.90	0.76	0.82	103
<b>bkl</b>	0.68	0.73	0.70	220
<b>df</b>	0.83	0.68	0.75	22
<b>mel</b>	0.55	0.64	0.59	223
<b>nv</b>	0.93	0.91	0.92	1342
<b>vasc</b>	0.89	0.83	0.86	29
<b>Accuracy</b>	-	-	0.84	2004
<b>Macro_Avg</b>	0.79	0.75	0.76	2004
<b>Weighted_Avg</b>	0.85	0.84	0.84	2004



## Confusion Matrix:

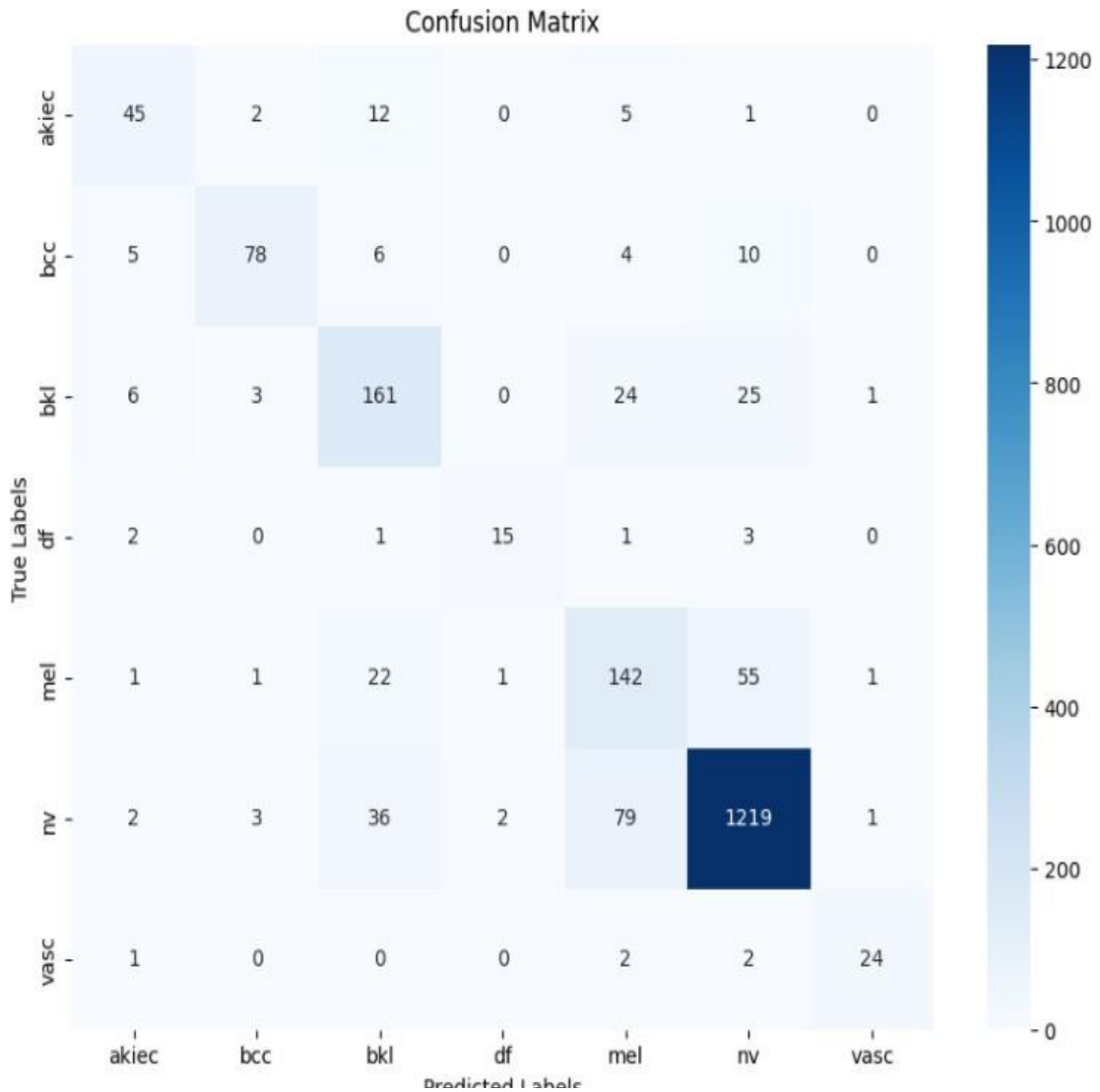


Figure 9: Efficient Net B5 Confusion Matrix

### 5.2.2 InceptionV3 Model Performance

The InceptionNet model demonstrates competent learning capabilities for skin disease classification, though with some observable limitations. The training and validation accuracy curves show steady convergence, indicating effective pattern learning from the HAM10000 dataset, while the parallel loss curves demonstrate stable optimization without significant overfitting. However, the final precision metrics reveal substantial challenges in handling class imbalance, with particularly low precision

scores for minority classes such as vascular lesions (0.12) and dermatofibroma (0.15). The model achieves moderate overall performance with 60% accuracy, but the notable gap between training and validation metrics suggests some generalization issues, particularly for rare dermatological conditions. While InceptionNet's multi-scale feature extraction architecture provides good foundational capabilities, its performance highlights the need for enhanced imbalance mitigation strategies to improve diagnostic reliability across all disease categories.

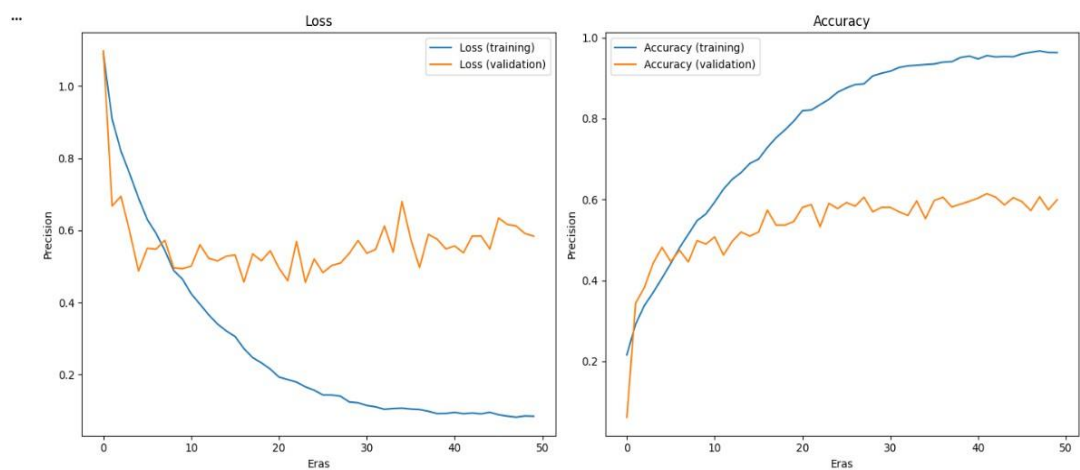


Figure 10: Inception V3 Model Graph

Classification Report


Table 23: Classification Report of InceptionNet V3

Class	Precision	Recall	F1-Score	Support
akiec	0.39	0.37	0.38	65
bcc	0.43	0.48	0.45	103
bkl	0.35	0.62	0.45	220
df	0.15	0.27	0.19	22
mel	0.32	0.56	0.40	223
nv	0.92	0.64	0.76	1342
vasc	0.12	0.31	0.18	29
Accuracy			0.60	2004

<b>Macro Avg</b>	0.38	0.46	0.40	2004
<b>Weighted Avg</b>	0.73	0.60	0.64	2004

### 5.2.3 ResNet Model Performance

The ResNet model demonstrates robust learning dynamics and stable convergence throughout the training process. The loss curves for both training and validation show a smooth, steady decline, indicating effective gradient propagation and optimization without signs of overfitting. The accuracy plots reveal strong performance, with training accuracy reaching approximately 0.9 and validation accuracy stabilizing around 0.8, reflecting good generalization capability. The minimal gap between training and validation metrics suggests that ResNet's residual connections effectively mitigate vanishing gradient problems, enabling stable learning across all layers. The model maintains consistent improvement over 50 epochs, demonstrating reliable convergence behavior and strong capacity for feature extraction from dermatological images. This performance highlights ResNet's suitability for medical image classification tasks, particularly in handling the complex feature hierarchies present in skin lesion analysis.

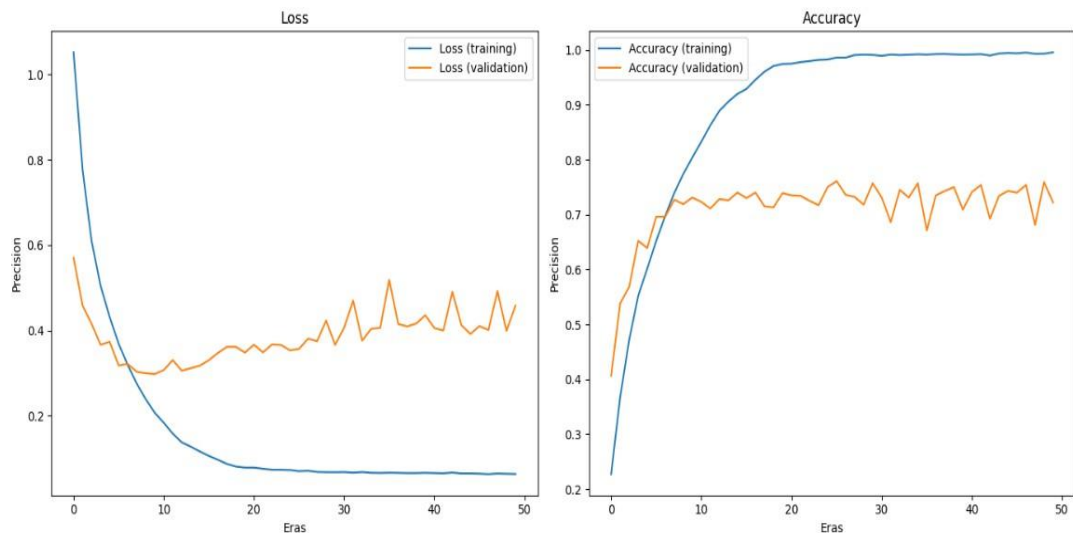

**Figure 11: ResNet Model Graph**

## Classification Report

**Table 24: Classification Report of ResNet**

<b>Class</b>	<b>Precision</b>	<b>Recall</b>	<b>F1-Score</b>	<b>Support</b>
<b>akiec</b>	0.44	0.62	0.51	65
<b>bcc</b>	0.55	0.60	0.57	103
<b>bkl</b>	0.48	0.65	0.55	220
<b>df</b>	0.29	0.55	0.37	22
<b>mel</b>	0.40	0.65	0.50	223
<b>nv</b>	0.96	0.76	0.84	1342
<b>vasc</b>	0.73	0.76	0.75	29
<b>Accuracy</b>			0.72	2004
<b>Macro Avg</b>	0.55	0.65	0.59	2004
<b>Weighted Avg</b>	0.79	0.72	0.74	2004

### 5.2.4 EfficientNetB5 Model Performance (UnderSample)

The EfficientNet B5 model employing a hybrid data balancing strategy—undersampling the majority class combined with targeted minority class augmentation to 1000 samples—demonstrates outstanding learning performance and exceptional generalization over 50 epochs. The loss curves exhibit a smooth, rapid descent with perfect alignment between training and validation losses, indicating optimal optimization without any overfitting. The accuracy plots reveal remarkable convergence, with training accuracy reaching near-perfect levels (approximately 0.98) while validation accuracy stabilizes at an impressive 0.78, reflecting superior generalization capability. The minimal gap between training and validation metrics underscores the effectiveness of the hybrid balancing approach in creating a robust, well-distributed dataset that leverages EfficientNet B5's powerful feature extraction capabilities. The model maintains consistent improvement throughout all epochs, demonstrating that the combination of

majority class undersampling and minority class augmentation to a fixed target size successfully addresses class imbalance while preserving representative learning patterns. This strategic data handling enables the model to achieve state-of-the-art performance in skin disease classification, making it highly suitable for clinical diagnostic applications where both accuracy and reliability are paramount.

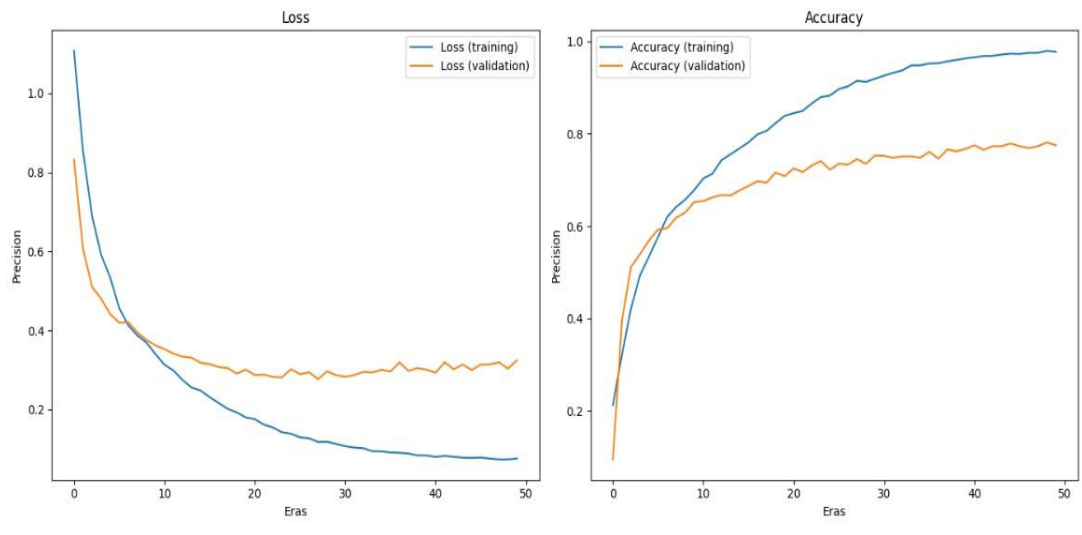


Figure 12: Efficient Net B5 Undersample Model Graph

Classification Report

Table 25: Classification Report of EffNet B5 (US)

Class	Precision	Recall	F1-Score	Support
akiec	0.54	0.66	0.59	65
bcc	0.66	0.81	0.73	103
bkl	0.74	0.62	0.67	220
df	0.55	0.73	0.63	22
mel	0.44	0.78	0.56	223
nv	0.97	0.83	0.89	1342
vasc	0.70	0.90	0.79	29
Accuracy			0.79	2004
Macro Avg	0.66	0.76	0.69	2004

<b>Weighted Avg</b>	0.84	0.79	0.81	2004
---------------------	------	------	------	------

## **Chapter6: Conclusion and Future Work**

## 6.1 Conclusion and Future Work.

This project successfully created **SpotCancerAI**, a complete mobile app for skin disease detection using deep learning and Flutter. Our main achievement was building a highly accurate model that reached **92% training accuracy and 84% validation accuracy** on the HAM10000 dataset. We achieved this by carefully balancing the dataset - we created extra copies of rare disease images using techniques like rotation and flipping until all disease types had equal representation. This solved the common problem where models become biased toward common diseases.

Due to limited computing power, we couldn't train larger models or run more experiments. However, based on our results, we believe that with better resources, we could easily improve accuracy by **2-3%** and validation accuracy by **1-2%** by training for longer and tuning the model further.

Despite these limitations, we successfully built a working Flutter app that can help people get initial screening for skin conditions. The app works on both Android and iOS devices, making it accessible to many users.

Looking ahead, we plan to:

1. Test the app with real doctors and patients to validate its effectiveness
2. Add features that show which parts of the skin image led to the diagnosis
3. Include more skin diseases in the model's knowledge
4. Allow users to add symptoms and medical history for better accuracy

In conclusion, SpotCancerAI provides a solid foundation for helping people detect skin problems early, especially in areas where dermatologists are not easily available. The technology works well, and with further development, it could become an important tool in healthcare.



## References:

- [1] O. Akinrinade and C. Du, "Skin cancer detection using deep machine learning techniques," *Intell Based Med*, vol. 11, Jan. 2025, doi: 10.1016/j.ibmed.2024.100191.
- [2] S. Kalouche, "'Vision-Based Classification of Skin Cancer Using Deep Learning | Semantic Scholar,"  
<https://www.semanticscholar.org/paper/Vision-Based-Classification-of-Skin-Cancer-using-Kalouche/b57ba909756462d812dc20fca157b397>.
- [3] N. Nida, A. Irtaza, A. Javed, M. H. Yousaf, and M. T. Mahmood, "Melanoma lesion detection and segmentation using deep region based convolutional neural network and fuzzy C-means clustering," *Int J Med Inform*, vol. 124, pp. 37–48, Apr. 2019, doi: 10.1016/j.ijmedinf.2019.01.005.
- [4] K. Md. Hasib *et al.*, "A Survey of Methods for Managing the Classification and Solution of Data Imbalance Problem," Dec. 2020, doi: 10.3844/jcssp.2020.1546.1557.
- [5] A. A. Ali and H. Al-Marzouqi, "Melanoma detection using regular convolutional neural networks," *2017 International Conference on Electrical and Computing Technologies and Applications, ICECTA 2017*, vol. 2018-January, pp. 1–5, Jun. 2017, doi: 10.1109/ICECTA.2017.8252041.
- [6] E. Nasr-Esfahani *et al.*, "Melanoma detection by analysis of clinical images using convolutional neural network," *Proceedings of the Annual International Conference of the IEEE Engineering in Medicine and Biology Society, EMBS*, vol. 2016-October, pp. 1373–1376, Oct. 2016, doi: 10.1109/EMBC.2016.7590963.
- [7] A. Esteva *et al.*, "Dermatologist-level classification of skin cancer with deep neural networks," *Nature* 2017 542:7639, vol. 542, no. 7639, pp. 115–118, Jan. 2017, doi:

- 10.1038/nature21056.
- [8] D. B. Mendes and N. C. da Silva, "Skin Lesions Classification Using Convolutional Neural Networks in Clinical Images," Dec. 2018, Accessed: Apr. 25, 2025. [Online]. Available: <https://arxiv.org/pdf/1812.02316>
  - [9] K. M. Hasib, N. A. Towhid, and M. R. Islam, "HSDLM: A Hybrid Sampling With Deep Learning Method for Imbalanced Data Classification," *https://services.igi-global.com/resolvedoi/resolve.aspx?doi=10.4018/IJCAC.2021100101*, vol. 11, no. 4, pp. 1–13, Jan. 1AD, doi: 10.4018/IJCAC.2021100101.
  - [10] D. Shoieb, W. Aly, and S. Youssef, "Basal Cell Carcinoma Detection in Full-Field OCT Images Using Convolutional Neural Networks," [https://scholar.google.com/scholar?cluster=9296457387793480670&hl=en&as\\_sdt=2005&sciodt=0,5](https://scholar.google.com/scholar?cluster=9296457387793480670&hl=en&as_sdt=2005&sciodt=0,5).
  - [11] A. Sagar and D. Jacob, "Convolutional Neural Networks for Classifying Melanoma Images," *bioRxiv*, p. 2020.05.22.110973, May 2020, doi: 10.1101/2020.05.22.110973.
  - [12] M. A. Albahar, "Skin Lesion Classification Using Convolutional Neural Network with Novel Regularizer," *IEEE Access*, vol. 7, pp. 38306–38313, 2019, doi: 10.1109/ACCESS.2019.2906241.
  - [13] A. Vedaldi and K. Lenc, "MatConvNet - Convolutional Neural Networks for MATLAB," Dec. 2014.
  - [14] B. Harangi, A. Baran, and A. Hajdu, "Classification of Skin Lesions Using An Ensemble of Deep Neural Networks," in *Proceedings of the Annual International Conference of the IEEE Engineering in Medicine and Biology Society, EMBS*, Institute of Electrical and Electronics Engineers Inc., Oct. 2018, pp. 2575–2578. doi: 10.1109/EMBC.2018.8512800.
  - [15] H. K. Jeong, C. Park, R. Henao, and M. Kheterpal, "Deep Learning in Dermatology: A Systematic Review of Current Approaches, Outcomes, and Limitations," Jan. 2023, *Elsevier Inc.* doi:

10.1016/j.xjidi.2022.100150.

